

Meicen Radiotherapy Positioning Solution



Quality · Service · Harmony · Innovation

Guangdong Meicen Medical Equipment Co., Ltd.

Tel: +86-757-86269115

E-mail: sales@meicenmed.com

Add: Building 5, Wanyang Innovation Park, No. 82 Xiaotang
Industrial Avenue, Shishan Town, Nanhai District, Foshan,
Guangdong, China.



Brochure 2024



www.meicenmed.com

Ver. 2402

Guangdong Meicen Medical Equipment Co., Ltd.



About Us



Guangdong Meicen Medical Equipment Co., Ltd. is an enterprise with independent R&D capabilities, and offering medical equipment R&D manufacturing, sales and services, supported with a more than 8000 square meters' modernized workshops and excellent R&D engineers, technicians, warehouse management and well-trained sales team.

Meicen is dedicated to providing high quality immobilization products for radiotherapy and orthopedic rehabilitation. Its core products include Thermoplastic Masks, Patient Positioning Devices and Components for Radiotherapy, Vacuum Cushions, Thermoplastic Splint Material, Thermoplastic Orthoses and Orthotic Components. All our products are in compliance with ISO13485:2016 certificates, SFDA, CE and other international regulations, and have been widely used and highly recognized, and hence Meicen has become one of the core suppliers of thermoplastic sheet materials in the world.

“Quality, Service, Harmony, Innovation” is our company philosophy, we take exquisite technique as core strength, keep improving and manufacturing supreme products with innovative technology, provide excellent service through scientific management together with our partners around the world.

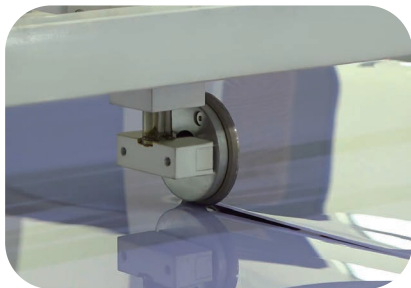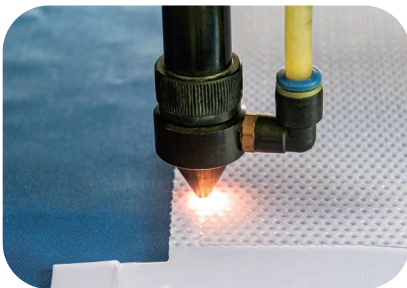

Contents

SBRT Immobilization System	2
C-Series Baseplates	8
Y-Series Baseplates	14
Supine Breast Cushions	18
A-Series Baseplates	20
B-Series Baseplates	23
Wing Board	24
Supine Breast Board	25
Prone Breast & Belly Solutions	26
Violet Thermoplastic Masks	28
Vacuum Cushions	50
Diagnostic Couchtop Overlays	53
Accessories	54

► SBRT Immobilization System

Meicen SBRT Immobilization System is highly accurate, comfortable, multi-functional and user-friendly device for SBRT treatment, which is compatible with the couchtops of Elekta, Varian or other brands. The paddle design on the belly and knee bridges can be adjusted in the direction of left and right, up and down to control the breathing range of the patient effectively during the treatment. The quick-lock and unlock system of the bridges is effortless. And, there are scales to record the fixed positions.

The system features a single-size breathing bridge that is fully adjustable to suit all patients. And the abundant scale marks improve the accuracy of positioning repeatability.

It can be used with Wing Board, Vacuum Cushions, Thermoplastic Masks, and other products to provide various immobilization solutions for Radiation Therapy.



SBRT Immobilization System Solutions

► SBRT Immobilization System



High accuracy

Clear and precise scale, improves the accuracy of positioning repeatability.

Patient comfort


Ergonomically design for patient comfort.

Versatility

Multifunctional design, simple to form a set-up for various demand.

User-friendly

Easy and simple to operate.

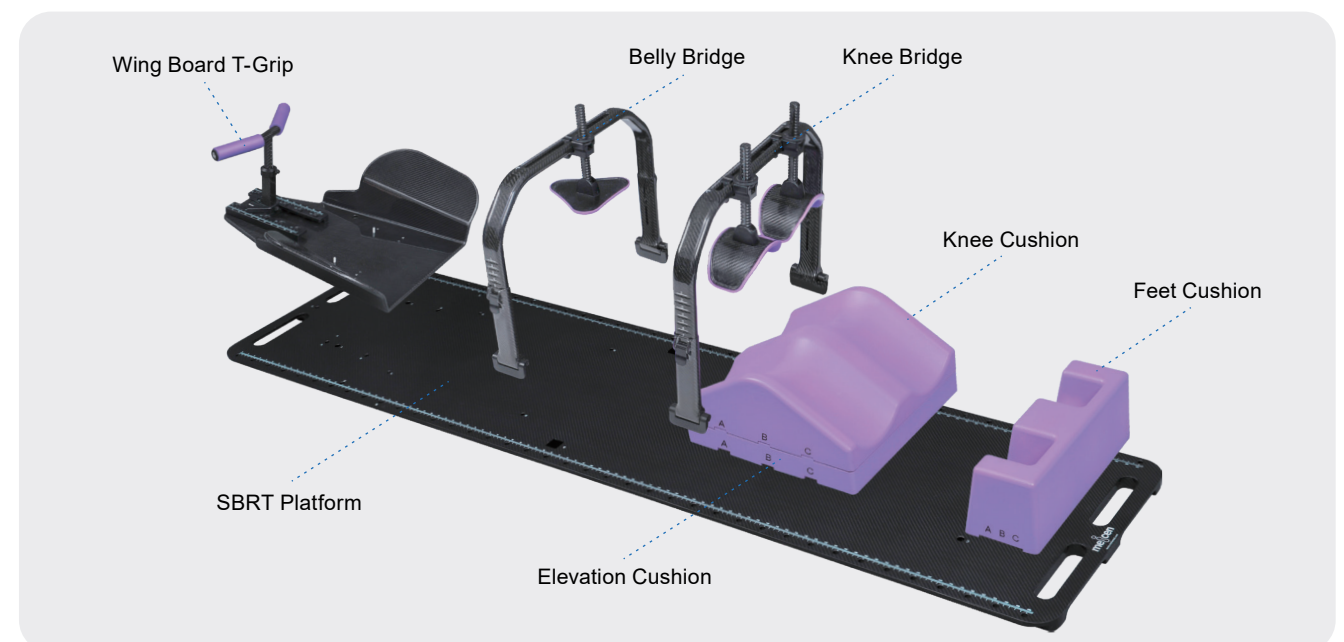

SBRT Immobilization System

REF MC-STF001




Carbon Fiber SBRT Immobilization System Package

	Components	REF	Qty
MC-STF001 Standard Package	SBRT Platform	MC-STF001-A01	1
	Wing Board T-Grip	MC-WBF001	1
	Belly Bridge	MC-STF001-A02	1
	Knee Cushion	RE-FCA01	1
	Elevation Cushion	RE-FCA02	1
	Feet Cushion	RE-FCA03	1
	T-Shaped Vacuum Cushion	RC-E3010F	1
	SBRT Indexing Bar	MC-STF001-A07	4
Optional Product	Carbon Fiber Headrest	RE-HSF6	1
	Color Headrest	RE-HSB5	1
	Clear Headrest	RE-HSC6	1
	Knee Bridge	MC-STF001-A03	1
	SRS Cradle	MC-SRF001	1
	Indexing Bar (for Elekta couchtop)	RC-LBE02	2
	Indexing Bar (for Varian couchtop)	RC-LBV02	2
	Shoulder Retractor	MC-STG001-A08	1

SBRT Immobilization System

REF MC-STG001

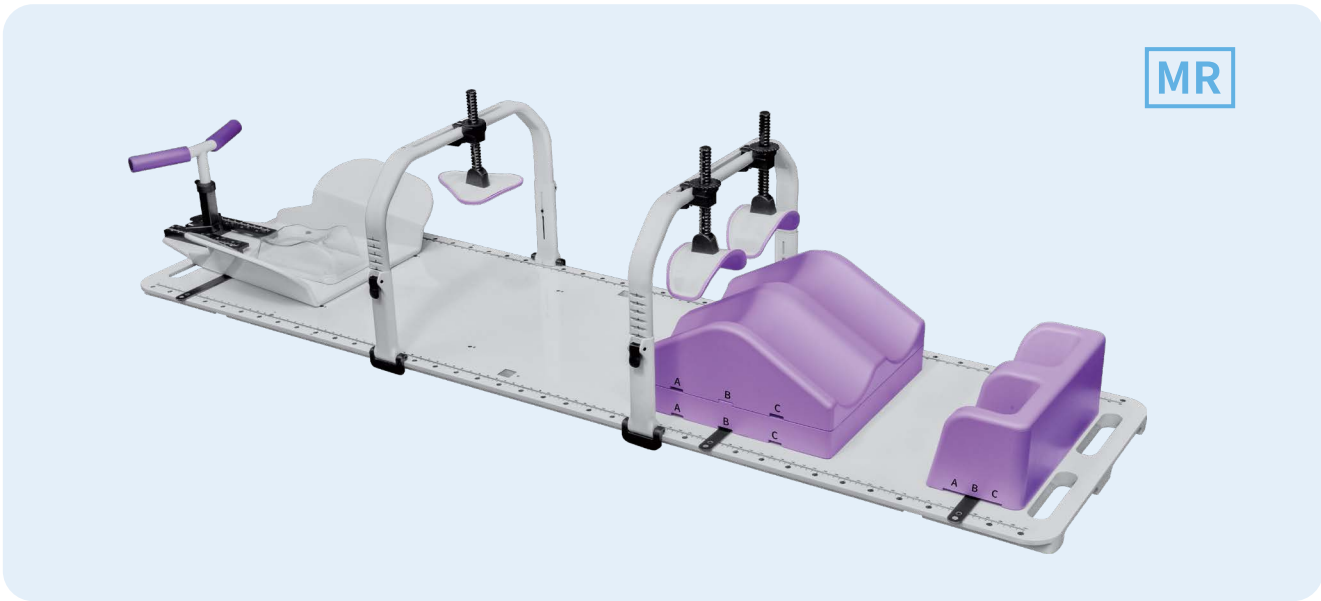

MR SBRT Immobilization System Package

	Components	REF	Qty
MC-STG001 Standard Package	SBRT Platform	MC-STG001-A01	1
	Wing Board T-Grip	MC-WBG001	1
	Belly Bridge	MC-STG001-A02	1
	Knee Cushion	RE-FCA01	1
	Elevation Cushion	RE-FCA02	1
	Feet Cushion	RE-FCA03	1
	T-Shaped Vacuum Cushion	RM-E3010F	1
	SBRT Indexing Bars	MC-STG001-A07	4
Optional Product	Color Headrest	RE-HSB5	1
	Clear Headrest	RE-HSC6	1
	Knee Bridge	MC-STG001-A03	1
	SRS Cradle	MC-SRG001	1
	Indexing Bar (for Elekta couchtop)	RM-LBE01	2
	Indexing Bar (for Varian couchtop)	RM-LBV01	2
	Shoulder Retractor	MC-STG001-A08	1


SBRT Immobilization System - Components

SBRT Platform

MC-STF001-A01 Carbon Fiber

MC-STG001-A01 Glass Fiber MR

- Specification: 2000mm (L) x 530mm (W).
- Multiple corresponding holes, efficient indexing, install without restraint.
- Lightweight, easy to move.

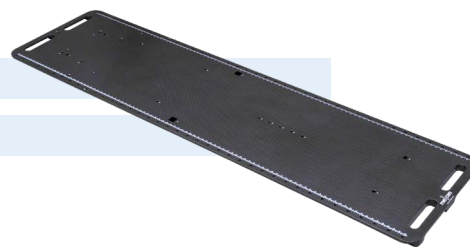

Wing Board T-Grip (Page 24)

MC-WBF001 Carbon Fiber

MC-WBG001 Glass Fiber MR

- T-Grip design provides comfortable arms extension.
- Can be used with the T-Shaped Vacuum Cushion, RC-E3005F or RC-E3010F.



Belly Bridge

MC-STF001-A02 Carbon Fiber

MC-STG001-A02 Glass Fiber MR

- Ergonomic Compression Paddle for better comfort.
- The Paddle device can be adjusted in the direction of left and right, up and down, and has scales to record its fixed position.
- Abundant scale marks improve the accuracy of positioning repeatability.



Knee Bridge

MC-STF001-A03 Carbon Fiber

MC-STG001-A03 Glass Fiber MR

- The knee paddle device can be adjusted left and right, up and down, and has scales to record its fixed position.
- Abundant scale marks improve the accuracy of positioning repeatability.
- Knee paddle device is optional in double leg or single leg immobilization.



SBRT Immobilization System - Components


Knee & Feet Cushion (Page 57)

RE-FCA01 Knee Cushion

RE-FCA03 Feet Cushion

RE-FCA02 Elevation Cushion

- Easy to operate, improve the comfort and stability of supine position during fixation.
- It can be directly fixed on the CT couchtop for use.

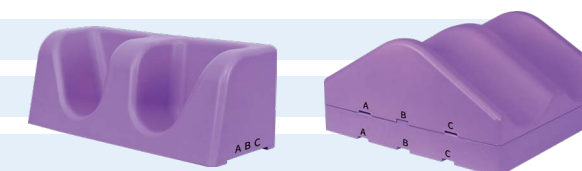

Feet Cushion

Knee & Elevation Cushion

T-Shaped Vacuum Cushion (Page 50-52)

RC-E3005F Non-open headrest

RC-E3010F Open headrest

- Specification: 125x110cm.
- Open headrest area can be used together with the headrest to reduce the ray blocking.
- Provides arms-up positioning and support the torso region.



Shoulder Retractor MC-STG001-A08

- It can be used with the SBRT Platform to improve patient comfort effectively and reduce involuntary displacement.



Headrest (Page 54)

RE-HSB5 Color Headrest, 5 Sizes, A to E

RE-HSC6 Clear Headrest, 6 Sizes, A to F

RE-HSF6 Carbon Fiber Headrest, 6 Sizes, A to F



RE-HSB5

RE-HSC6


RE-HSF6

SBRT Indexing Bar

MC-STF001-A07 Carbon Fiber

MC-STG001-A07 Glass Fiber MR


- Allows positioning devices to be fixed to the SBRT Platform for enhancing reproducibility.

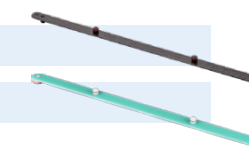

Meicen C-Series Baseplates

The Meicen C-Series Baseplates have a Pin-lock design for easy operation, provides the highest degree of stability for head, head-shoulders, thorax, pelvic and full-body support. Dual-configured for adult and pediatric treatments. The C-Series Baseplates are lightweight, and adjustable shoulder retractor with shoulder suppression. It is compatible with all S-type Masks and Meicen Pin-lock type Masks.

- Contoured design for patients comfort
- Sides emargination for simple installation and carry
- Configured for multiple treatments
- Different location choices to provide better immobilization for each patient

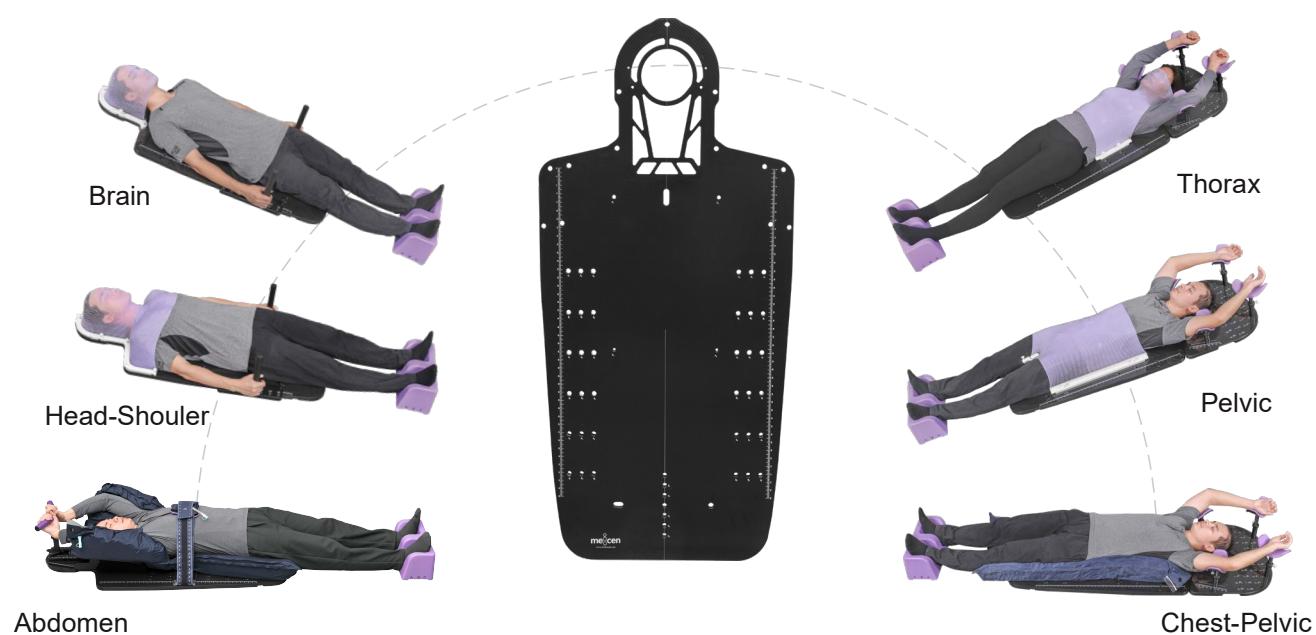


C-Series AIO Baseplate Solutions

Meicen C-Series Baseplates


Meicen C-Series AIO Baseplate

MC-CF001	Carbon Fiber
MC-CG001	Glass Fiber MR


- Provide multiple treatment options for head, head-shoulders, thorax, pelvic and full-body support.
- Different location choices to provide better immobilization for each patient.
- The Leg Separator can be locked and removed easily, and adjusted to fit each patient.
- It can be adapted to the standard and narrow shoulder S-type Masks.
- Optional Overhead Arm Support and Shoulder Retractor can be used to improve the accuracy and comfort during treatment.
- Specification: 1300(L) x 606(W)/mm

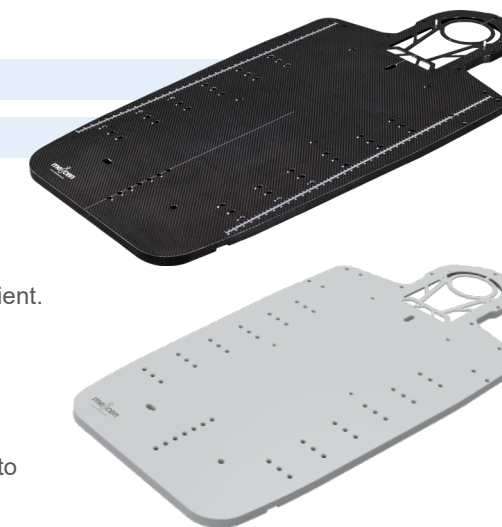

Easily adjustable Hand Grips

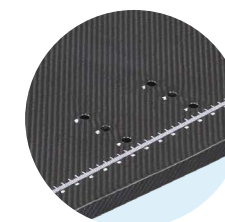



Adjustable to patient size

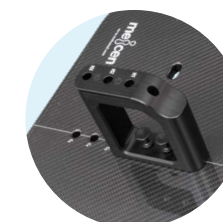

Leg Separator (optional)

Extension Board RC-EBF01

- Light weight, easy to move.
- Offers advanced indexing options.
- MR compatible: RM-EBG01

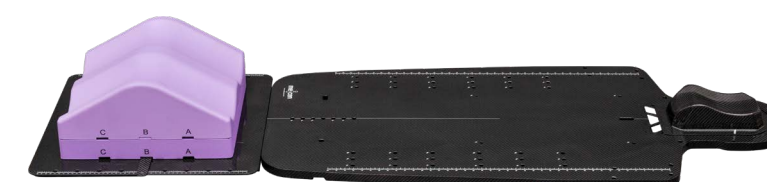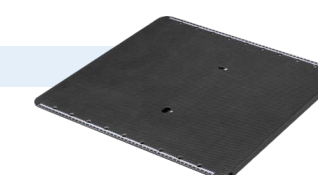

Meicen C-Series Baseplates

Meicen C-Series Head-Shoulder-Breast Baseplate

MC-CF002 Carbon Fiber

MC-CG002 Glass Fiber MR

- The contoured design for patients comfort.
- Sides emargination for simple installation and carry.
- Provide the highest degree of stability for head, head-shoulders and breast.
- It can be adapted to the standard and narrow shoulder S-type Masks.
- Comfortable and adjustable Hand Grips assist with shoulder retraction.
- Specification: 1000(L) x 606(W)/mm




Meicen C-Series Head-Shoulder Baseplate

MC-CF004 Carbon Fiber

MC-CG004 Glass Fiber MR

- Specific design for Elekta HexaPOD™ evo RT System.
- Dual-configured for standard and narrow shoulder S-type Masks.
- Specification: 550(L) x 550(W)/mm

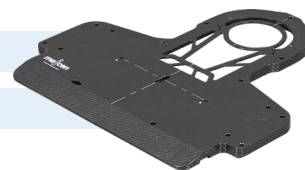

Meicen C-Series Pelvic Baseplate

MC-CF003 Carbon Fiber

MC-CG003 Glass Fiber MR

- The contoured design for patients comfort.
- Sides emargination for simple installation and carry.
- Different location choices to provide better immobilization for each patient.
- Specification: 850(L) x 575(W)/mm



Meicen C-Series Baseplates - Optional Components

Shoulder Retractor can effectively improve patient comfort and reduce inactive displacement, comfortable and adjustable Hand Grips assist with shoulder retraction.

Pin-Lock Shoulder Retractor MC-CG001-A01

Compatible with MC-CF001 and MC-CG001



Pin-Lock Shoulder Retractor MC-CG002-A02

Compatible with MC-CF002 and MC-CG002



Can be locked and removed easily, and adjusted to fit each patient.

Leg Separator MC-CG001-A04

Compatible with MC-CF001 and MC-CG001



Groin Lock MC-CG001-A03



Compatible with MC-CF001, MC-CG001, MC-CF003 and MC-CG003



Meicen C-Series Baseplates - Optional Components

S-Shaped Head Mask (Page 32&34)

- M2422-VSA190W 22% perforation, 2.4mm
- M3038-VSA480W 38% perforation, 3.0mm
- M24RA-VSA580W Reinforced, 2.4mm
- M30RA-VSA610W Reinforced, 3.0mm

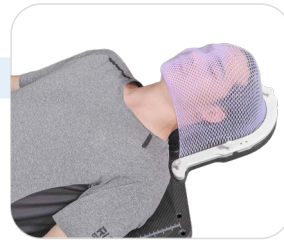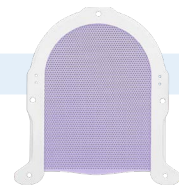


S-Shaped Head-Neck-Shoulder Mask (Page 33&35)

- M2422-VSB190W 22% perforation, 2.4mm
- M3038-VSB480W 38% perforation, 3.0mm
- M24RA-VSB580W Reinforced, 2.4mm
- M30RA-VSB610W Reinforced, 3.0mm
- M2422-VSB191W 22% perforation, 2.4mm, narrow shoulder




Pin-lock type Neck-Thorax Mask (Page 42)

- M24RA-VSF580W Reinforced, 2.4mm
- M30RA-VSF610W Reinforced, 3.0mm

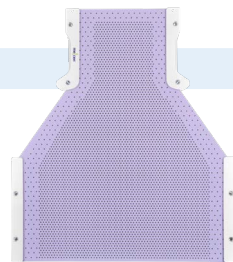

Pin-lock type Trunk Mask (Chest-Pelvic) (Page 42)

- M2422-VST190W 22% perforation, 2.4mm
- M3038-VST480W 38% perforation, 3.0mm

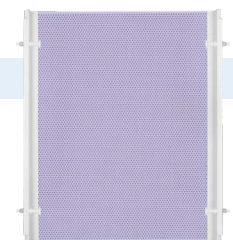

Pin-lock type Trunk Mask (Pelvic) (Page 42)

- M2422-VST191W 22% perforation, 2.4mm
- M3038-VST481W 38% perforation, 3.0mm



Meicen C-Series Baseplates - Optional Components

Vacuum Cushion (Page 50-52)

- 100x70cm, with Indexing Bars, RC-E5005S
- 150x75cm, with Indexing Bars, RC-E5033S
- 100x70cm, RC-E5003
- 150x75cm, RC-E5032



S-Type Headrests (Page 54)

- RE-HSF6 Carbon Fiber Headrest, 6 Sizes, A to F
- RE-HSB5 Color Headrest, 5 Sizes, A to E
- RE-HSC6 Clear Headrest, 6 Sizes, A to F



RE-HSF6




RE-HSB5



RE-HSC6

Knee & Feet Cushion

- Improve the comfort of the patient when positioning.
- Can be fixed directly to the CT or treatment couchtop, or with the Meicen Extension Board.

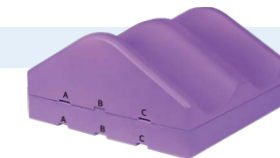

Knee Cushion RE-FCA01



Feet Cushion RE-FCA03

Elevation Cushion RE-FCA02

Overhead Arm Support RC-ASF01

- It can be used together with the Meicen C-Series Baseplate MC-CF001 and MC-CF002 to improve reproducible precision and comfort in thorax and lung tumor treatment.
- The height of Wrist Supports can be adjusted to each patient.
- The angle of Arms Supports can be adjusted to each patient.
- MR compatible: RM-ASG01



Meicen Y-Series Baseplates

The Meicen Y-Series Baseplates have a Pin-lock design for easy operation, provide the highest degree of stability for head, head-shoulders, thorax, pelvic and full-body support. Dual-configured for adult and pediatric treatments. It is compatible with all S-type Masks and Meicen Pin-lock type Masks. The Y-Series AIO Baseplate can be combined with Meicen Belly & Knee Bridge for SBRT treatment needs.

- Contoured design for patients comfort
- Compatible with narrow shoulder S-type Masks
- Configured for multiple treatments



Meicen Y-Series Baseplates

Meicen Y-Series AIO Baseplate

MC-YF001	Carbon Fiber
MC-YG001	Glass Fiber MR

- Provide multiple treatment options for head, head-shoulders, thorax, pelvic and full-body support.
- Adaptable to the standard and narrow shoulder S-type Masks.
- Optional Overhead Arm Support and Shoulder Retractor can be used to improve the accuracy and comfort during treatment.
- Specification: 1300(L) x 606(W)/mm




Extended SBRT Solution for Y-Series AIO Baseplate

Components	CT Package	MR Package	Qty
Belly Bridge	MC-YF001-A01	MC-YG001-A01	1
Knee Bridge	MC-YF001-A02	MC-YG001-A02	1
Extension Board	RC-EBF01	RM-EBG01	1
Overhead Arm Support	RC-ASF01	RM-ASG01	1
Knee Cushion	RE-FCA01	RE-FCA01	1
Elevation Cushion	RE-FCA02	RE-FCA02	1
Feet Cushion	RE-FCA03	RE-FCA03	1
Headrest	RE-HSF6	RE-HSB5	1
Indexing Bar (for Elekta couchtop)	RC-LBE02	RM-LBE01	2
Indexing Bar (for Varian couchtop)	RC-LBV02	RM-LBV01	2
Shoulder Retractor	MC-YG001-A03	MC-YG001-A03	1

Meicen Y-Series Baseplates - Optional Components

S-Shaped Head Mask (Page 32&34)

- M2422-VSA190W 22% perforation, 2.4mm
- M3038-VSA480W 38% perforation, 3.0mm
- M24RA-VSA580W Reinforced, 2.4mm
- M30RA-VSA610W Reinforced, 3.0mm

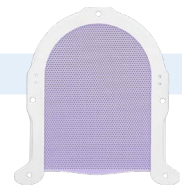


S-Shaped Head-Neck-Shoulder Mask (Page 33&35)

- M2422-VSB190W 22% perforation, 2.4mm
- M3038-VSB480W 38% perforation, 3.0mm
- M24RA-VSB580W Reinforced, 2.4mm
- M30RA-VSB610W Reinforced, 3.0mm
- M2422-VSB191W 22% perforation, 2.4mm, narrow shoulder




Pin-lock type Neck-Thorax Mask (Page 42)

- M24RA-VSF580W Reinforced, 2.4mm
- M30RA-VSF610W Reinforced, 3.0mm

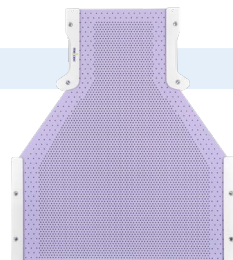

Pin-lock type Trunk Mask (Chest-Pelvic) (Page 42)

- M2422-VST190W 22% perforation, 2.4mm
- M3038-VST480W 38% perforation, 3.0mm

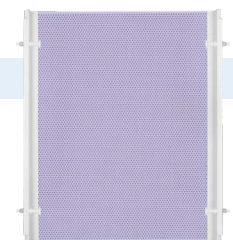

Pin-lock type Trunk Mask (Pelvic) (Page 42)

- M2422-VST191W 22% perforation, 2.4mm
- M3038-VST481W 38% perforation, 3.0mm



Meicen Y-Series Baseplates - Optional Components

Vaccum Cushion (Page 50-52)

- 100x80cm with Indexing Bar, RC-E5009F
- 150x80cm with Indexing Bars, RC-E5026S
- 100x80cm, RC-E5016
- 150x80cm, RC-E5027




Overhead Arm Support RC-ASF01

- It can be used together with the MC-YF001 to improve the reproducible precision and comfort in breast thorax and lung tumor treatment.
- The height of Wrist Supports and angle of Arm Supports can be adjusted to each patient.
- MR compatible: RM-ASG01



Extension Board RC-EBF01

- Improve the comfort of the patient when positioning.
- Can be fixed directly to the CT or treatment couchtop.
- MR compatible: RM-EBG01

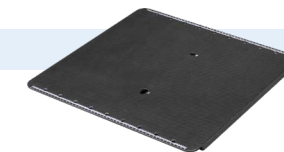

Pin-Lock Shoulder Retractor MC-YG001-A03

- Compatible with MC-YF001 and MC-YG001



Belly Bridge MC-YF001-A01


- Carbon fiber material.
- Ergonomic Compression Paddle for better comfort.
- The Paddle device can be adjusted in the direction of left and right, up and down, and has a scale to record its fixed position.
- Abundant scale marks improve the accuracy of positioning repeatability.




Supine Breast Cushions

Supine Breast Cushions MC-RXE001

The Supine Breast Cushions can be easily installed on Meicen C & Y Series Baseplates by adapter cushion, or can also be indexed directly to couchtops. The Arm Support Cushions is ergonomically designed for patient comfort. It can be used alone or in conjunction with the Meicen C & Y Series Baseplates. The Hand Grips can be adjusted easily to improve the accuracy and comfort of repeatability.

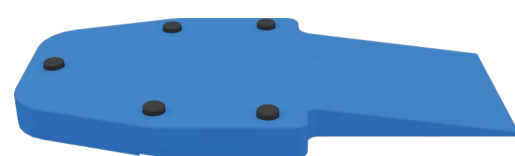


Wedge Cushion 8°
MC-RXE001-A01

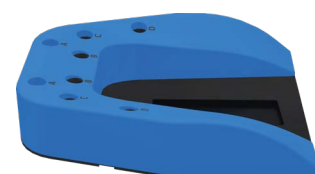

Arm Support Cushion (A-type)
MC-RXE001-A03



Wedge Cushion 13°
MC-RXE001-A02

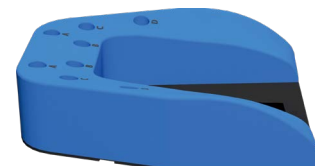

Arm Support Cushion (B-type)
MC-RXE001-A04



Supine Breast Cushions



Ergonomic design for patient comfort



Indexed slot for fixing on CT couchtops




Quick-lock design to easily lock into
Meicen baseplates



The better positioning locations according
to patient condition

Freedom Solution

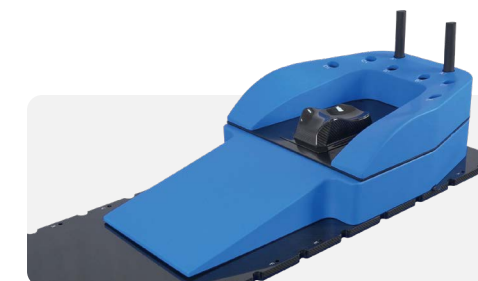

Please contact our sales representative for more solutions.


Meicen A-Series Baseplates

Meicen A-Series AIO Baseplate

MC-AF001 Carbon Fiber

MC-AG001 Glass Fiber MR

- Provides multiple treatment options for head, neck, shoulders, thorax, pelvic and full-body support.
- Different location choices to provide better immobilization for each patient.
- Shoulder Retractor can be used to improve the accuracy and comfort during treatment.
- Specification: 1300(L) x 550(W)/mm

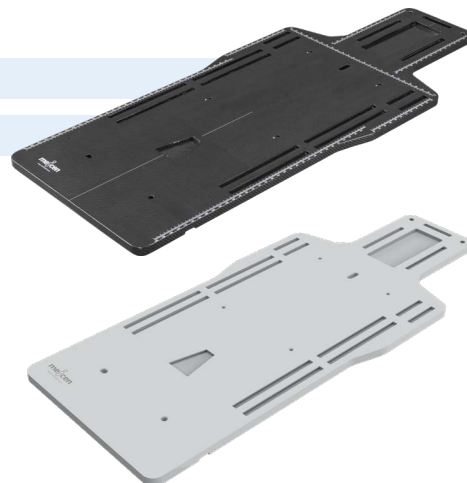

Shoulder Retractor MC-AG001-A01

- Comfort and reduce inactive displacement effectively.
- Compatible with MC-AF001 and MC-AG001 for better patient comfort.



Compatible with



Vacuum Cushions - page 50 to 52

Meicen A-Series Baseplates

Meicen A-Series Head-Shoulder-Breast Baseplate

MC-AF002 Carbon Fiber

MC-AG002 Glass Fiber MR

- Contoured design for patient comfort.
- Sides emargination for simple installation and carry.
- Provides immobilization solution for head, neck, shoulders and breast.
- Specification: 1000(L) x 550(W)/mm



Meicen A-Series Pelvic Baseplate

MC-AF003 Carbon Fiber

MC-AG003 Glass Fiber MR

- Contoured design for patients comfort.
- Sides emargination for simple installation and carry.
- Provides immobilization solution for thorax, abdomen and pelvic.
- Specification: 850(L) x 520(W)/mm




Meicen A-Series Head-Shoulder Baseplate

MC-AF004 Carbon Fiber

MC-AG004 Glass Fiber MR


- Provides immobilization solution for head, neck and shoulders.
- Compatible with all U-frame Masks and L-type Masks.
- Compatible with Color Headrest RE-HSP5.
- Specification: 550(L) x 550(W)/mm

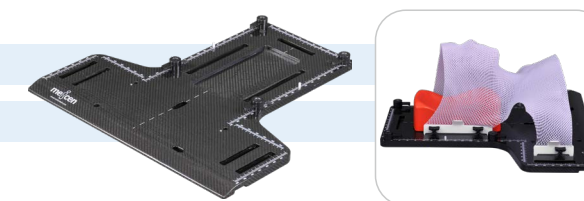

Meicen A-Series Baseplates - Optional Components


Triangular Head Mask (Page 36)

- M2422-VLA191W 22% perforation, 2.4mm
- M24RA-VLA580W Reinforced, 2.4mm

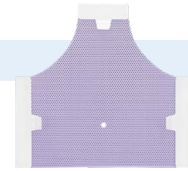


Head-Neck-Shoulder Mask (Page 37)

- M2422-VLB191W 22% perforation, 2.4mm
- M24RA-VLB580W Reinforced, 2.4mm

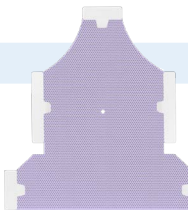

Chest-Pelvic Mask (Page 37)

- M3038-VLK481W 38% perforation, 3.0mm

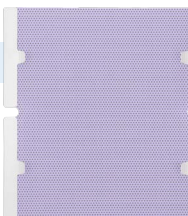


Vacuum Cushion (Page 50-52)

- 120x80cm, with Indexing Bars, RC-E5025S
- 100x75cm, RC-E5020
- 70x70cm, RC-E4013



Color Headrest RE-HSP5


- Effectively improve patient comfort and reduce involuntary displacement.
- 5 Sizes, A to E.
- Can be compatible with A-Series Baseplates.

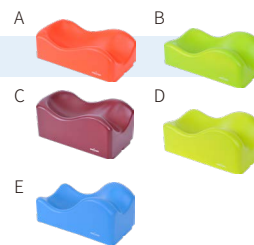

Meicen A-Series Baseplates - Optional Components

Block Set (standard)

- MC-AF001-A02 Carbon Fiber
- MC-AG001-A03 MR

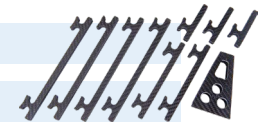

Hand Grips (standard)

- RC-HGF01 Carbon Fiber
- RM-HGG01 MR




Meicen B-Series Baseplates

Meicen B-Series U-Frame Head Baseplate

- MC-BF003 Carbon Fiber
- MC-BG003 Glass Fiber MR


- Provides the head support.
- Compatible with all U-frame Masks.
- Compatible with Meicen all type headrests.
- Specification: 366(L) x 420(W)/mm

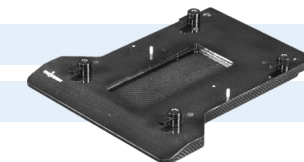

Meicen B-Series Tilting U-Frame Head Baseplate

- MC-BF004 Carbon Fiber
- MC-BG004 Glass Fiber MR

- Provides the head support.
- Compatible with all U-frame Masks.
- Tilting angles: 5°/10°/20°/25°/30°/35°/40°/45°.
- Compatible with Meicen all type headrests.
- Specification: 377(L) x 306(W)/mm

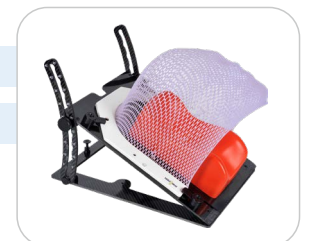

Meicen Wing Board

MC-WBF001	T-Grip, Carbon Fiber
MC-WBF002	U-Grip, Carbon Fiber

- The Meicen Wing Board is integrally formed with high-quality carbon fiber material, which is durable and lightweight.
- It's easy to use with Vacuum Cushion greatly improves the comfort of the patient's treatment and the accuracy of positioning.
- The patient's arm is raised to maintain a grasping state, which can improve the stability and accuracy of the patient's immobilization.
- The Meicen Wing Board can be used alone or in conjunction with the Meicen C-Series Baseplates, SBRT System and Supine Breast Board.



Wing Board T-Grip



Wing Board U-Grip

MC-WBF001-A01	T-Grip Carbon Fiber
MC-WBG001-A01	T-Grip MR



RC-HGF01	U-Grip Carbon Fiber
RM-HGG01	U-Grip MR



MC-WBG001 (MR Series)



T-Shaped Vacuum Cushion RC-E3010F

Meicen Supine Breast Board

MC-RXF002	Carbon Fiber
MC-RXG002	Glass Fiber MR

- The Supine Breast Board has 5 inclinations: 5°, 10°, 15°, 20° and 25°.
- Simple, customizable setups for patient comfort and clearance.
- Contoured lumbar area and enhanced bottomstop with easy locking.
- Compatible with S-type Head Masks for stable, reproducible supraclavicular treatments.
- Compatible with Meicen Breast Masks.
- Specification: 1325(L) x 445(W)/mm



Arm Support MC-RXF001-A02




Wrist Support MC-RXF001-A01



Compatible with

M2422-VSH190W	Breast Mask, 22% perforation, 2.4mm
M3038-VSH480W	Breast Mask, 38% perforation, 3.0mm

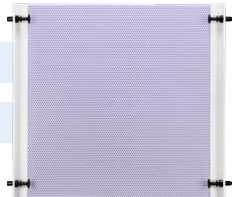

Compatible with Breast Mask can be reduced the patient's involuntary displacement and improve the repeatability and comfort of the positioning.



Meicen Prone Breast & Belly Solutions

Meicen Prone Breast Board

MC-RBF001	Carbon Fiber
MC-RBG001	Glass Fiber MR

- Variable aperture openings to accommodate different breast sizes.
- Multiple heterolateral wedge sizes to assist with patient tilt.
- Adjustable and durable massage-style face cushion.
- Scale rulers for easy verification of longitudinal and lateral position of the patient.
- Reversible breast section allows for treatment of either breast.
- Specification: 1570 (L) x 450 (W) x 150 (H)/mm



MC-RBG001-A01



MC-RBG001-A02



MC-RBG001-A03



MC-RBG001-A04




MC-RBG001-A05

Compatible with

M2422-VHT289W	Trunk Mask, 22% perforation, 2.4mm
M3038-VHT550W	Trunk Mask, 38% perforation, 3.0mm

Compatible with Trunk Mask together can be reduced the patient's involuntary displacement and improve the repeatability.

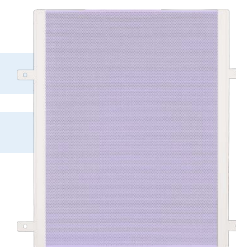

Meicen Prone Breast & Belly Solutions

Prone Breast Cushion MC-RBE001

- The Breast-bridge can be converted left and right according to the treatment plan.
- The U-shaped cushion can be selected S, M, L size.
- The rail-type prone headrest with head to feet direction and the four-levels elevation angle can be adjusted to meet individual needs.
- Compatible with C & Y-Series AIO Baseplate.



Prone Headrest



Prone Belly Cushion MC-RPE002

- The rail-type prone headrest with head to feet direction and the four-levels elevation angle can be adjusted to meet individual needs.
- Dual-configured for female and male treatments.
- It can be used with Trunk Mask to improve the immobilization effect.
- Compatible with C & Y-Series AIO Baseplate.



Prone Headrest



Meicen Violet Thermoplastic Masks


Meicen Violet Thermoplastic Masks are composed of a new, premium thermoplastic material designed specifically to offer the ultimate in patient comfort with exceptional rigidity and extremely low shrinkage. Available in multiple perforation styles and compatible with most manufacturers' baseplates, our violet thermoplastic's primary focus is the patient experience.

- Low shrinkage
- Optimal molding time
- Ergonomically designed for patient comfort
- Non-Stick
- Violet color to alleviate patient anxiety




Meicen Violet Thermoplastic Masks


Meicen Reinforced Series Masks feature a combination of enhanced perforation and solid thermoplastic to provide rigidity and secure immobilization for IMRT treatments. Openface Masks can be compatible with SGRT (Surface Guided Radiation Therapy) camera systems.

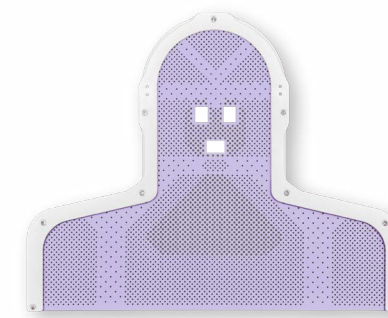


S-type Reinforced

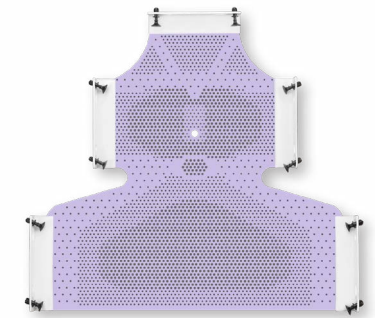

P-type Reinforced

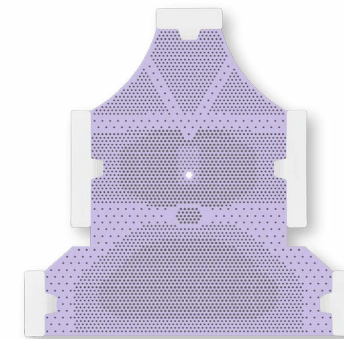

L-type Reinforced

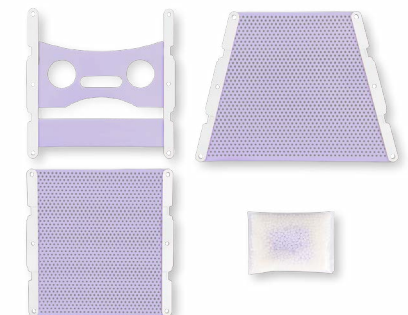

Stereotactic X-Knife

Meicen Violet Thermoplastic Masks are easy to prepare and can be heated by water bath at 65-70°C (149-158°F) or oven at 70-75°C (158-167°F).

Moulding Time: 1-2 minutes

Cooling Time: 8-15 minutes

Thickness 2.0mm, 2.4mm, 3.0mm or 3.2mm are available.

Please contact our sales representative for more information.

► U-frame Violet Thermoplastic Masks

- Slim Frame



U-Shaped Head Mask

Standard


M2422-VUA192W	22% perforation, 2.4mm
M3038-VUA481W	38% perforation, 3.0mm
M2415-VUA150W	15% perforation, 2.4mm



U-Shaped Head Mask

Reinforced


M24RA-VUA580W	2.4mm
M30RA-VUA610W	3.0mm

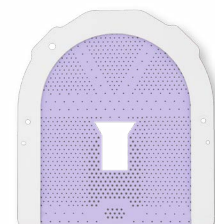

U-Shaped Head Mask

Reinforced, Openface

M24RB-VUA640W	2.4mm
M30RB-VUA670W	3.0mm

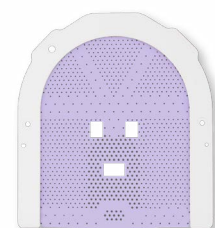

U-Shaped Head Mask

Reinforced, Open eyes & mouth

M24RB-VUA641W	2.4mm
M30RB-VUA671W	3.0mm



U-Shaped Head Mask

Standard, Extended

Extended 5cm	Extended 9cm	
M2422-VUA197W	M2422-VUA200W	22% perforation, 2.4mm
M3038-VUA482W	M3038-VUA483W	38% perforation, 3.0mm

► U-frame Violet Thermoplastic Masks

- Standard Frame



U-Shaped Head Mask

Standard


M2422-VUA190W	22% perforation, 2.4mm
M3038-VUA480W	38% perforation, 3.0mm



U-Shaped Head Mask

Reinforced


M24RA-VUA581W	2.4mm
M30RA-VUA611W	3.0mm

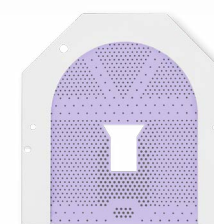

U-Shaped Head Mask

Reinforced, Openface


M24RB-VUA642W	2.4mm
M30RB-VUA672W	3.0mm

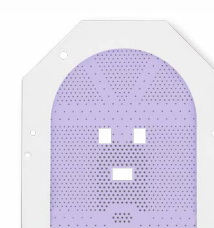

U-Shaped Head Mask

Reinforced, Open eyes & mouth

M24RB-VUA643W	2.4mm
M30RB-VUA673W	3.0mm

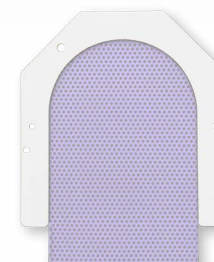

U-Shaped Head Mask

Standard, Extended

Extended 5cm	Extended 9cm	
M2422-VUA194W	M2422-VUA202W	22% perforation, 2.4mm
M3038-VUA484W	M3038-VUA485W	38% perforation, 3.0mm

Meicen U-frame Masks fit all standard U-frame baseplates from Meicen and other manufacturers.

2.0mm & 3.2mm thickness are available also, please specify thickness when requesting a quote.

► S-type Violet Thermoplastic Masks



S-Shaped Head Mask

Standard

M2422-VSA190W	22% perforation, 2.4mm
M3038-VSA480W	38% perforation, 3.0mm
M2415-VSA150W	15% perforation, 2.4mm



S-Shaped Head Mask

Reinforced

M24RA-VSA580W	2.4mm
M30RA-VSA610W	3.0mm



S-Shaped Head Mask

Standard, Extended

Extended 5cm		Extended 9cm
M2422-VSA192W	M2422-VSA193W	22% perforation, 2.4mm
M3038-VSA482W	M3038-VSA483W	38% perforation, 3.0mm



S-Shaped Head Mask

Reinforced, Extended

Extended 5cm		Extended 9cm
M24RA-VSA581W	M24RA-VSA582W	2.4mm
M30RA-VSA611W	M30RA-VSA612W	3.0mm

Meicen S-type Masks fit all standard S-type baseplates from Meicen and other manufacturers.
2.0mm & 3.2mm thickness are available also, please specify thickness when requesting a quote.


► S-type Violet Thermoplastic Masks



S-Shaped Head-Neck-Shoulder Mask

Standard


M2422-VSB190W	22% perforation, 2.4mm
M3038-VSB480W	38% perforation, 3.0mm
M2415-VSB150W	15% perforation, 2.4mm

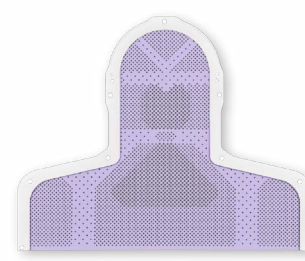

S-Shaped Head-Neck-Shoulder Mask

Reinforced


M24RA-VSB580W	2.4mm
M30RA-VSB610W	3.0mm

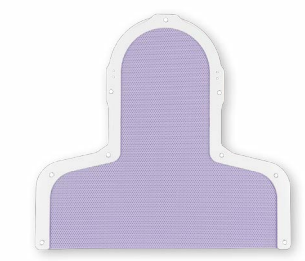

S-Shaped Head-Neck-Shoulder Mask (narrow shoulders)

Standard

M2422-VSB191W	22% perforation, 2.4mm
M3038-VSB481W	38% perforation, 3.0mm

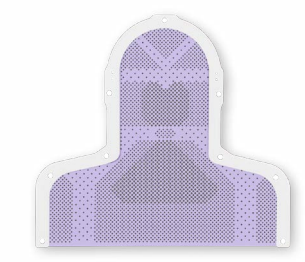

S-Shaped Head-Neck-Shoulder Mask (narrow shoulders)

Reinforced


M24RA-VSB581W	2.4mm
M30RA-VSB611W	3.0mm

* Compatible with Meicen C & Y-Series baseplates.

► S-type Violet Thermoplastic Masks (Openface)




Meicen Openface Series Masks aims to alleviate symptoms of claustrophobia during radiotherapy. All types of masks are available as openfaced.
Adjustable pulling of the open part when apply to different patients assures comfort and reinforced effect.

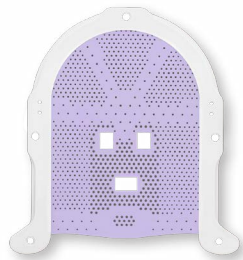


S-Shaped Head Mask
Reinforced, Open eyes & mouth

M24RB-VSA643W	2.4mm
M30RB-VSA672W	3.0mm

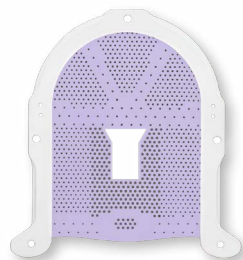

S-Shaped Head Mask
Reinforced, Openface

M24RB-VSA640W	2.4mm
M30RB-VSA670W	3.0mm

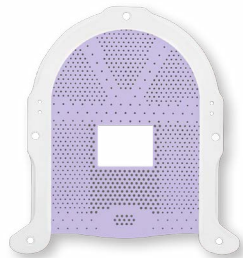

S-Shaped Head Mask
Reinforced, Openface

M24RB-VSA641W	2.4mm
M30RB-VSA671W	3.0mm

► S-type Violet Thermoplastic Masks (Openface)



S-Shaped Head-Neck-Shoulder Mask
Reinforced, Openface

M24RB-VSB642W	2.4mm
M30RB-VSB670W	3.0mm



S-Shaped Head-Neck-Shoulder Mask
Reinforced, Open eyes & mouth

M24RB-VSB643W	2.4mm
M30RB-VSB672W	3.0mm



S-Shaped Head-Neck-Shoulder Mask
Reinforced, Openface

M24RB-VSB641W	2.4mm
M30RB-VSB671W	3.0mm



S-Shaped Head-Neck-Shoulder Mask
Reinforced, Open face & neck

M24RB-VSB644W	2.4mm
M30RB-VSB673W	3.0mm

Meicen S-type Masks fit all standard S-type baseplates from Meicen and other manufacturers.
2.0mm & 3.2mm thickness are available also, please specify thickness when requesting a quote.

► L-type Violet Thermoplastic Masks



Triangular Head Mask

Standard, with nose hole

M2422-VLA191W	22% perforation, 2.4mm
M3038-VLA480W	38% perforation, 3.0mm
M2438-VLA320W	38% perforation, 2.4mm, for pediatric



Triangular Head Mask

Standard, Extended 3cm, with nose hole


M2422-VLA192W	22% perforation, 2.4mm
M3038-VLA481W	38% perforation, 3.0mm
M2415-VLA151W	15% perforation, 2.4mm



Triangular Head Mask

Reinforced, with nose hole


M24RA-VLA580W	2.4mm
M30RA-VLA610W	3.0mm

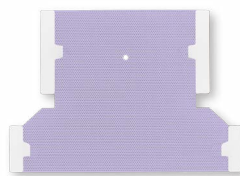

Neck-Shoulder Mask

Standard, with nose hole

M2422-VLE190W	22% perforation, 2.4mm
M3038-VLE480W	38% perforation, 3.0mm
M2422-VLE191W	22% perforation, 2.4mm, for pediatric

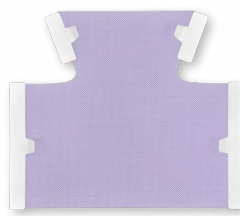


Neck-Thorax Mask

Standard

M2422-VLF192W	22% perforation, 2.4mm
M3038-VLF480W	38% perforation, 3.0mm


► L-type Violet Thermoplastic Masks

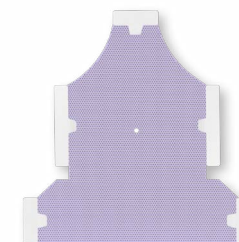

Head-Neck-Shoulder Mask

Standard, with nose hole


M2422-VLB191W	22% perforation, 2.4mm
M3038-VLB481W	38% perforation, 3.0mm
M2422-VLB192W	22% perforation, 2.4mm, for pediatric

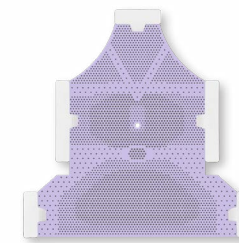

Head-Neck-Shoulder Mask

Reinforced, with nose hole


M24RA-VLB580W	2.4mm
M30RA-VLB610W	3.0mm

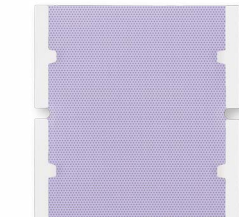

Chest-Pelvic Mask

Standard


M3038-VLK483W	38% perforation, 3.0mm, 228 clamp
M3038-VLK482W	38% perforation, 3.0mm, 243 clamp
M3038-VLK481W	38% perforation, 3.0mm, 247 clamp

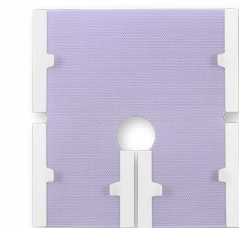

Pelvic Mask

Standard

M3038-VLL483W	38% perforation, 3.0mm, 228 clamp
M3038-VLL482W	38% perforation, 3.0mm, 243 clamp
M3038-VLL481W	38% perforation, 3.0mm, 247 clamp

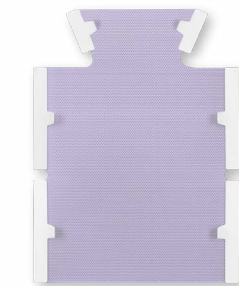


Thorax-Abdomen Mask

Standard

M2422-VLG190W	22% perforation, 2.4mm
M3038-VLG480W	38% perforation, 3.0mm

Meicen L-type Masks fit all standard O-type baseplates from Meicen and other manufacturers.
2.0mm & 3.2mm thickness are available also, please specify thickness when requesting a quote.


► L-type Violet Thermoplastic Masks (Openface)

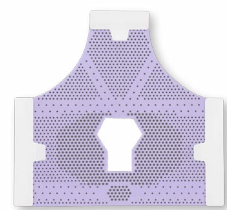

Triangular Head Mask

Reinforced, Openface


M24RB-VLA641W	2.4mm
M30RB-VLA671W	3.0mm

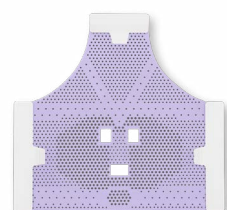

Triangular Head Mask

Reinforced, Open eyes & mouth

M24RB-VLA643W	2.4mm
M30RB-VLA673W	3.0mm

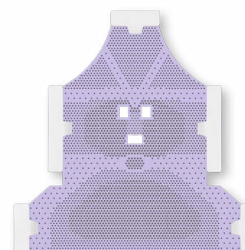

Head-Neck-Shoulder Mask

Reinforced, Open eyes & mouth


M24RB-VLB642W	2.4mm
M30RB-VLB672W	3.0mm



Head-Neck-Shoulder Mask

Openface

M24RB-VLB643W	2.4mm
M30RB-VLB673W	3.0mm

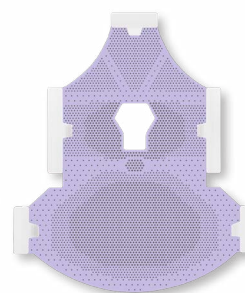

Head-Neck-Shoulder Mask

Reinforced, Extended, Openface

M24RB-VLC641W	2.4mm
M30RB-VLC671W	3.0mm

Meicen L-type Masks fit all standard O-type baseplates from Meicen and other manufacturers.
2.0mm & 3.2mm thickness are available also, please specify thickness when requesting a quote.

► P-type Violet Thermoplastic Masks



Triangular Head Mask

Standard, with nose hole


M2422-VPA190W	22% perforation, 2.4mm
M3038-VPA480W	38% perforation, 3.0mm
M2415-VPA150W	15% perforation, 2.4mm



Triangular Head Mask

Reinforced, with nose hole


M24RA-VPA580W	2.4mm
M30RA-VPA610W	3.0mm

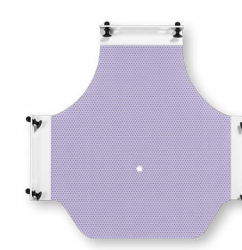

Triangular Head Mask

Standard, Extended 5cm, with nose hole


M2422-VPA191W	22% perforation, 2.4mm
M3038-VPA481W	38% perforation, 3.0mm
M2415-VPA151W	15% perforation, 2.4mm

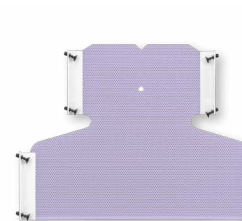

Neck-Shoulder Mask

Standard, with nose hole

M2422-VPE190W	22% perforation, 2.4mm
M3038-VPE480W	38% perforation, 3.0mm

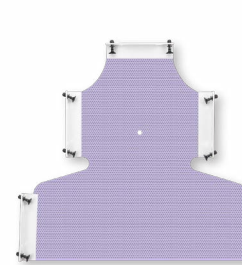


Head-Neck-Shoulder Mask

Standard, with nose hole

M2422-VPB190W	22% perforation, 2.4mm
M3038-VPB480W	38% perforation, 3.0mm
M2415-VPB150W	15% perforation, 2.4mm

► P-type Violet Thermoplastic Masks

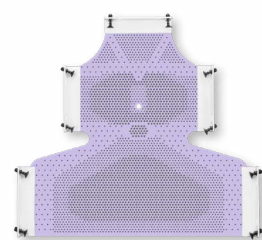

Head-Neck-Shoulder Mask

Reinforced, with nose hole

M24RA-VPB580W	2.4mm
M30RA-VPB610W	3.0mm



Chest-Pelvic Mask

Standard

M2422-VPK190W	22% perforation, 2.4mm
M3038-VPK481W	38% perforation, 3.0mm



Triangular Head Mask (for pediatric)

Standard, with nose hole


M2422-VPA192W	22% perforation, 2.4mm
M3038-VPA482W	38% perforation, 3.0mm



Head-Neck-Shoulder Mask (for pediatric)

Standard, with nose hole

M2422-VPB191W	22% perforation, 2.4mm
M3038-VPB481W	38% perforation, 3.0mm

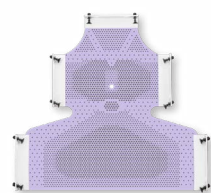

Head-Neck-Shoulder Mask (for pediatric)

Reinforced, with nose hole

M24RA-VPB581W	2.4mm
M30RA-VPB611W	3.0mm

Meicen P-type Masks fit all PosiFix-type baseplates from Meicen and Civco.

2.0mm & 3.2mm thickness are available also, please specify thickness when requesting a quote.

► P-type Violet Thermoplastic Masks (Openface)



Triangular Head Mask

Standard, Openface

M24RB-VPA640W	2.4mm
M30RB-VPA670W	3.0mm



Triangular Head Mask

Reinforced, Openface


M24RB-VPA641W	2.4mm
M30RB-VPA671W	3.0mm



Triangular Head Mask

Reinforced, Open eyes & mouth


M24RB-VPA642W	2.4mm
M30RB-VPA672W	3.0mm

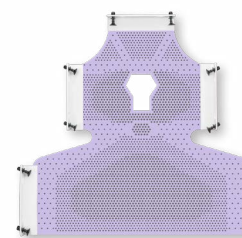

Head-Neck-Shoulder Mask

Reinforced, Openface

M24RB-VPB640W	2.4mm
M30RB-VPB670W	3.0mm

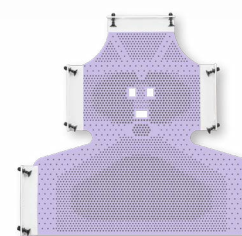


Head-Neck-Shoulder Mask

Reinforced, Open eyes & mouth

M24RB-VPB641W	2.4mm
M30RB-VPB671W	3.0mm

Pin-lock type Violet Thermoplastic Masks

- For the Meicen C & Y-Series baseplates.

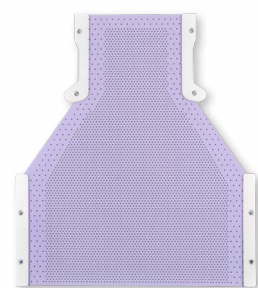

Neck-Thorax Mask

Reinforced

M24RA-VSF580W	2.4mm
M30RA-VSF610W	3.0mm



Trunk Mask (Chest-Pelvic)

Standard

M2422-VST190W	22% perforation, 2.4mm
M3038-VST480W	38% perforation, 3.0mm




Trunk Mask (Pelvic)

Standard, with Groin Lock

M2422-VST191W	22% perforation, 2.4mm
M3038-VST481W	38% perforation, 3.0mm

For Meicen Supine Breast Board

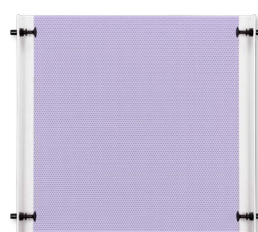

Breast Mask


Standard

M2422-VSH190W	22% perforation, 2.4mm
M3038-VSH480W	38% perforation, 3.0mm

2.0mm & 3.2mm thickness are available also, please specify thickness when requesting a quote.

H-type Violet Thermoplastic Masks


- For the Civco HipFix System.

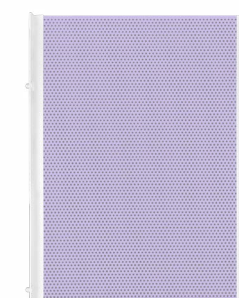

Trunk Mask (Chest-Pelvic)

Standard


M2422-VHT190W	22% perforation, 2.4mm
M3038-VHT480W	38% perforation, 3.0mm

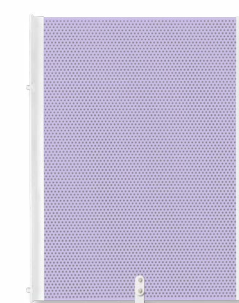

Trunk Mask (Pelvic)

Standard, with Groin Lock


M2422-VHT191W	22% perforation, 2.4mm
M3038-VHT481W	38% perforation, 3.0mm

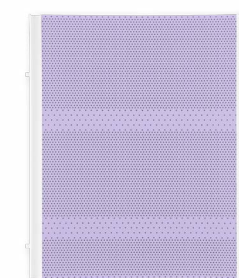

Trunk Mask (Chest-Pelvic)

Reinforced

M24RA-VHT580W	2.4mm
M30RA-VHT610W	3.0mm

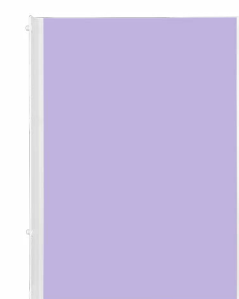

Trunk Mask (Chest-Pelvic)


Plain

M2400-VHT120W	2.4mm
M3000-VHT420W	3.0mm

2.0mm & 3.2mm thickness are available also, please specify thickness when requesting a quote.

M-type Violet Thermoplastic Masks


- For the MacroMedics OmniBoard 2.0

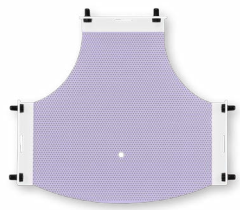

Triangular Head Mask

Standard, with nose hole

M2422-VMA190W	22% perforation, 2.4mm
M3038-VMA480W	38% perforation, 3.0mm

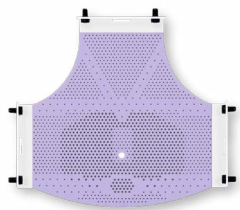

Triangular Head Mask

Reinforced, with nose hole


M24RA-VMA580W	2.4mm
M30RA-VMA610W	3.0mm



Head-Neck-Shoulder Mask

Standard, with nose hole


M2422-VMB190W	22% perforation, 2.4mm
M3038-VMB480W	38% perforation, 3.0mm

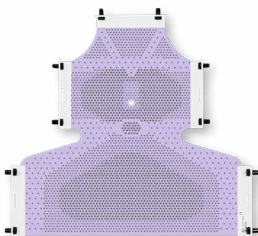

Head-Neck-Shoulder Mask

Reinforced, with nose hole

M24RA-VMB580W	2.4mm
M30RA-VMB610W	3.0mm

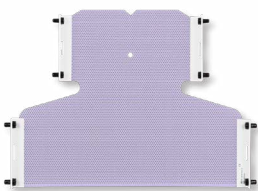

Neck-Shoulder Mask


Standard

M2422-VME190W	22% perforation, 2.4mm
M3038-VME480W	38% perforation, 3.0mm

X-Knife Violet Thermoplastic Masks

Meicen X-Knife Violet Thermoplastic Masks are composed of a new, premium thermoplastic material designed specifically to offer the ultimate in patient comfort with exceptional rigidity and extremely low shrinkage. Additionally, the X-Knife Violet Thermoplastic Masks have a non-stick coating applied to the patient-facing mask surfaces, helping ensure that the masks do not stick to the patient's hair, beard, etc.

- Adhesive Layers
- Low shrinkage
- Easy to manipulate and mould
- Cost-effective

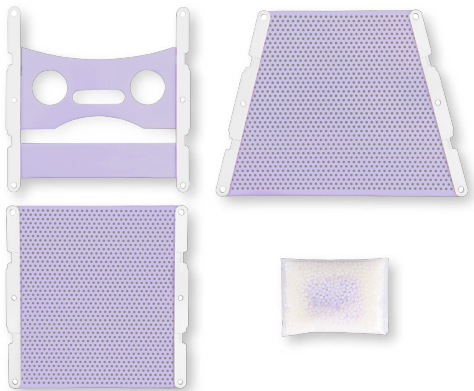


X-Knife Mask

38% perforation, 3.0mm
M3038-VXA481W



X-Knife Head-Shoulder Mask

38% perforation, 3.0mm
M3038-VXB480W

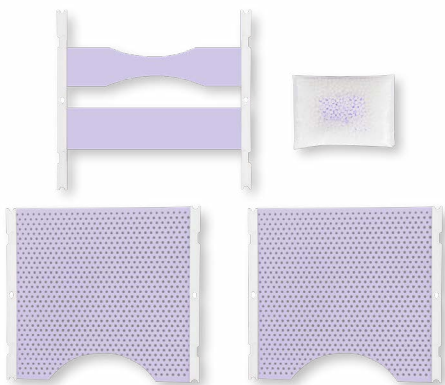

X-Knife Mask

38% perforation, 3.0mm
M3038-VXA490W

2.0mm & 3.2mm thickness are available also, please specify thickness when requesting a quote.

► E-type Violet Thermoplastic Masks

Compatible with Elekta Fraxion System



Triangular Head Mask

Reinforced

M24RA-VEA600W	2.4mm
M30RA-VEA620W	3.0mm



Triangular Head Mask

Reinforced

M24RA-VEA602W	2.4mm
M30RA-VEA621W	3.0mm

Compatible with Elekta HeadSTEP & HipSTEP Baseplates



Head Mask

Standard

M2422-VEA190W	22% perforation, 2.4mm
M3038-VEA480W	38% perforation, 3.0mm
M2415-VEA150W	15% perforation, 2.4mm




Head Mask

Reinforced

M24RA-VEA580W	2.4mm
M30RA-VEA610W	3.0mm

2.0mm & 3.2mm thickness are available also, please specify thickness when requesting a quote.

► E-type Violet Thermoplastic Masks

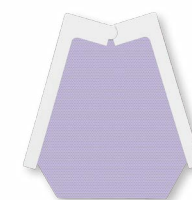

Head Mask

Standard, Extended 5cm


M2422-VEA192W	22% perforation, 2.4mm
M3038-VEA482W	38% perforation, 3.0mm



Head Mask

Standard, Extended 21cm


M2422-VEA193W	22% perforation, 2.4mm
M3038-VEA483W	38% perforation, 3.0mm

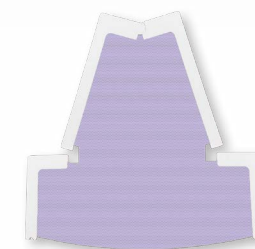

Head-Neck-Shoulder Mask

Standard

M2422-VEB190W	22% perforation, 2.4mm
M3038-VEB480W	38% perforation, 3.0mm
M2415-VEB150W	15% perforation, 2.4mm

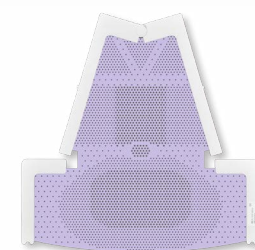

Head-Neck-Shoulder Mask

Reinforced

M24RA-VEB580W	2.4mm
M30RA-VEB610W	3.0mm



Chest-Pelvic Mask

Standard

M3038-VEK480W	38% perforation, 3.0mm
---------------	------------------------

► V-type Violet Thermoplastic Masks



V-Shaped Head Mask

Standard


M2422-VVA190W	22% perforation, 2.4mm
M3038-VVA480W	38% perforation, 3.0mm



V-Shaped Head Mask

Reinforced

M24RA-VVA580W	2.4mm
M30RA-VVA610W	3.0mm

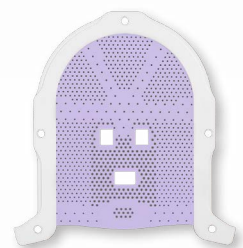

V-Shaped Head Mask

Reinforced, Open eyes & mouth

M24RB-VVA640W	2.4mm
M30RB-VVA670W	3.0mm



V-Shaped Head Mask

Reinforced, Openface

M24RB-VVA641W	2.4mm
M30RB-VVA671W	3.0mm



V-Shaped Head Mask

Standard, Extended

Extended 5cm		Extended 9cm
M2422-VVA191W	M2422-VVA192W	22% perforation, 2.4mm
M3038-VVA481W	M3038-VVA482W	38% perforation, 3.0mm

► V-type Violet Thermoplastic Masks



V-Shaped Head-Neck-Shoulder Mask

Standard


M2422-VVB190W	22% perforation, 2.4mm
M3038-VVB480W	38% perforation, 3.0mm



V-Shaped Head-Neck-Shoulder Mask

Reinforced

M24RA-VVB580W	2.4mm
M30RA-VVB610W	3.0mm

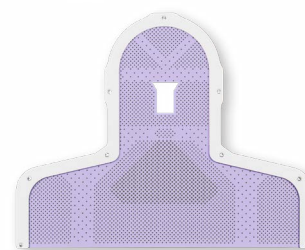

V-Shaped Head-Neck-Shoulder Mask

Reinforced, Openface

M24RB-VVB640W	2.4mm
M30RB-VVB670W	3.0mm



V-Shaped Head-Neck-Shoulder Mask

Reinforced, Open eyes & mouth

M24RB-VVB641W	2.4mm
M30RB-VVB671W	3.0mm

For the Klarity Optek board and Bionix VersaBoard.

2.0mm & 3.2mm thickness are available also, please specify thickness when requesting a quote.

Meicen Vacuum Cushions

The Meicen Vacuum Cushions offer a rigid and secure support around the patients. The coated nylon shell protects them from tear and stain. The Vacuum Cushions are completely air-tight, leak proof, easy clean and fill with optional density particles made from polyamide. Built-in indexing bars are available, different size and shape of Vacuum Cushions are available too.



Meicen Vacuum Cushions

- Ultra durable, high-quality polyamide.
- Dark / Light blue color distinguishable in optical imagers.
- Radio-translucent.
- Easy to clean.
- Cost-effective.
- Adapters available for any manufacturer's pump.



Valve



Indexing Bar



CPC Adapter

CT	MR	Dimensions	Fill	Indexed	
RC-E5003 RC-E5005S	RM-E5003 RM-E5005S	100x70cm	30L	W/Indexed	For Chest-Pelvic Support
RC-E4013 RC-E4014F	RM-E4013 RM-E4014F	70x70cm	20L	W/Indexed	For Breast Support
RC-E5024 RC-E5025S	RM-E5024 RM-E5025S	120x80cm	32L	W/Indexed	For Chest-Pelvic Support
RC-E5010 RC-E5012S	RM-E5010 RM-E5012S	100x75cm	30L	W/Indexed	For Chest-Pelvic Support
RC-E5023S	RM-E5023S	120x75cm	34L	W/Indexed	For Chest-Pelvic Support
RC-E5032 RC-E5033S	RM-E5032 RM-E5033S	150x75cm	45L	W/Indexed	For Chest-Pelvic Support

Remark: Custom width shall not exceed 140cm
Custom length shall not exceed 250cm
please specify color (Dark blue or Light blue color) when requesting a quote.

Meicen Vacuum Cushions



T-Shape Vacuum Cushion
with Open Headrest, 125 x 110cm
RC-E3010F



T-Shape Vacuum Cushion
125 x 110cm
RC-E3005F



Wing Board Cushion
110 x 55cm
RC-E4003



S-Shape Vacuum Cushion
61 x 65cm
RC-E1004

Non-Rectangular Vacuum Cushions

CT	MR	Dimensions	Fill	Shape	
RC-E1005	RM-E1005	75x75cm	15L	Arms up	For Breast Support
RC-E3005F RC-E3010F	RM-E3005F RM-E3010F	125x110cm	31L 32L	T-Shaped	For Thorax-Abdomen Support
RC-E3019	RM-E3019	105x75cm	28L	Arms up	For Thorax-Abdomen Support
RC-E4003 RC-E4004	RM-E4003 RM-E4004	110x55cm	15L 19L	Wing board breast	For Breast Support
RC-E1004	RM-E1004	61x65cm	8L	S-Shaped	For Head-Shoulder Support

Diagnostic Couchtop Overlays

- The Couchtop Overlays provide a clean, flat surface for accurate and easily reproducible patient positioning.
- The internal density distribution of the product is uniform, and there is no unacceptable artifact in the scanning area.
- Can be quickly installed and disassembled, simple and convenient.
- It is stable and reliable after installation and meets the leveling requirements.
- The Couchtop Overlays are designed to work with Indexing bars as well as other positioning systems.


* The Couchtop Overlays can be customized according to the user's CT or MRI equipment.

MC-CTF019 CT Overlay



Compatible to Somatom CT Definition AS+

MC-MRG002 MRI Overlay

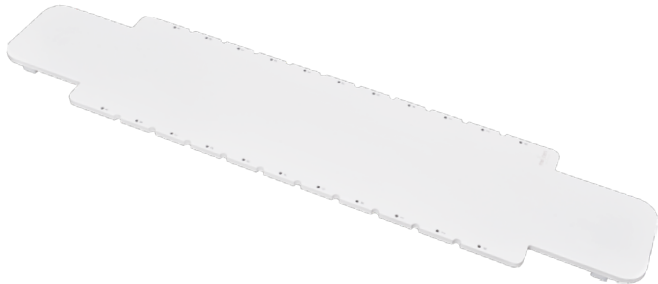


Compatible to Philips Ingenia Overlay

Headrests

The headrests included in sets can also be ordered separately, and for other manufacturers' baseplates are available upon request.

Carbon Fiber Headrest RE-HSF6

- Sturdy and durable.
- 6 sizes A to F.
- Compatible with all S-type baseplates.

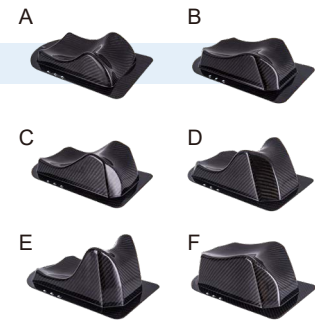


Color Headrest, S-type RE-HSB5

- PU material with customized CT or MR safe, soft and comfortable.
- 5 sizes A to E.
- Compatible with all S-type baseplates.




Clear Headrest RE-HSC6

- 6 sizes A to F.
- Compatible with all S-type baseplates.

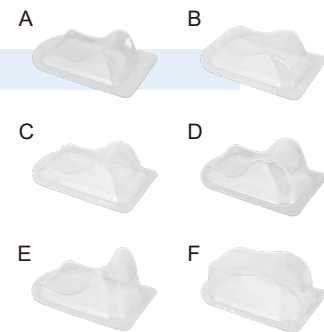

Color Headrest, A-type RE-HSP5

- Customized CT or MR safe, soft and comfortable.
- 5 sizes A to E.
- Compatible with Meicen A & P Series and MacroMedics Baseplates.

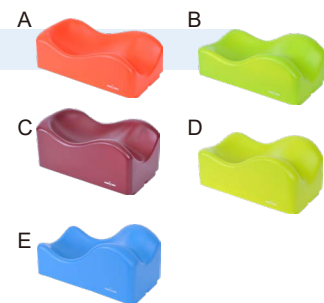

Shim


- Shim for Head, RE-HSA04
- Shim for Head-Neck-Shoulder, RE-HSA05



Indexing Bars


For Elekta Couchtops

- RC-LBE01 Aluminum
- RC-LBE02 Carbon Fiber
- RM-LBE01 Glass Fiber MR

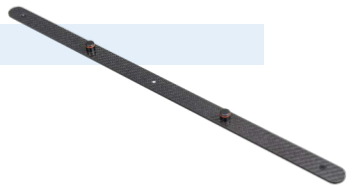


For Varian Couchtops

- RC-LBV01 Aluminum
- RC-LBV02 Carbon Fiber
- RM-LBV01 Glass Fiber MR

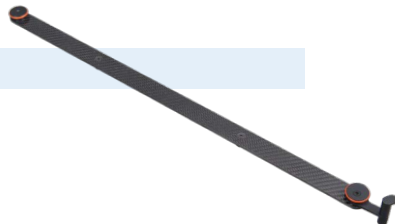

For Vacuum Cushions

- RC-LBC07 Carbon Fiber

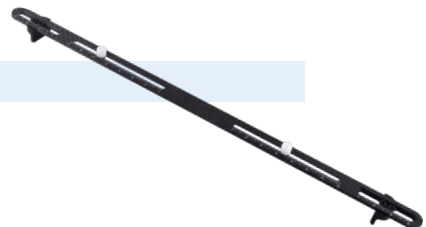

For Meicen C-Series AIO Baseplate


- RC-LBV03 Carbon Fiber for Varian Couchtops
- RC-LBE03 Carbon Fiber for Elekta Couchtops



Vacuum Cushion Accessories

Vacuum Cushion Positioners Strip RC-LBC02

- Carbon Fiber

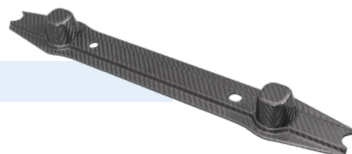

Vacuum Pump R-PM002

- 220V voltage
- It has two functions: vacuuming and inflating with two mechanical pressure gauges.



Others

Pin-Lock Shoulder Retractor MC-CG001-A01

- It can be used with the MC-CF001 and MC-CG001 which can effectively improve patient comfort and reduce inactive displacement.



Pin-Lock Shoulder Retractor MC-CG002-A02

- Compatible with MC-CF002 and MC-CG002.



Pin-Lock Shoulder Retractor MC-YG001-A03

- Compatible with Meicen MC-YF001 and MC-YG001.



Shoulder Retractor MC-AG001-A01

- It can be used together with MC-AF001 and MC-AG001 to effectively improve patient comfort and reduce inactive displacement.



Hand Grips RC-HGF01 Carbon Fiber

Hand Grips RM-HGG01 MR

- It can be used with the Meicen Shoulder Retractors and Baseplates.




Knee Cushion RE-FCA01

- Easy to operate, improve the comfort and stability of supine position during fixation.
- It can be directly fixed on the CT or MR couchtop for leg immobilization.



Elevation Cushion RE-FCA02

- It can be used together with the Knee Cushion RE-FCA01.

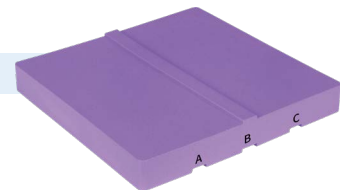

Feet Cushion RE-FCA03

- Easy to operate, improve the comfort and stability of supine position during fixation.
- It can be directly fixed on the CT or MR couchtop for feet immobilization.



Water Bath 04RF-G2

- Inner size: 700 x 550 x 120mm.
- Stainless steel material, suitable for 220V voltage.



Bolus

- 300 x 300 x 5mm
- 300 x 300 x 10mm



Purchasing Index

SBRT Immobilization System

REF	Product Description	Pkg	Page
MC-STF001	SBRT Immobilization System Carbon Fiber	SET	2-7
MC-STG001	SBRT Immobilization System Glass Fiber MR	SET	2-7

Meicen Y-Series Baseplates

REF	Product Description	Pkg	Page
MC-YF001	Y-Series AIO Baseplate Carbon Fiber	EA.	15
MC-YG001	Y-Series AIO Baseplate Glass Fiber MR	EA.	15
MC-YF001-A01	Belly Bridge Carbon Fiber	EA.	17
MC-YG001-A01	Belly Bridge Glass Fiber MR	EA.	17
MC-YF001-A02	Knee Bridge Carbon Fiber	EA.	17
MC-YG001-A02	Knee Bridge Glass Fiber MR	EA.	17
MC-YG001-A03	Pin-Lock Shoulder Retractor	EA.	17

Meicen C-Series Baseplates

REF	Product Description	Pkg	Page
MC-CF001	C-Series AIO Baseplate Carbon Fiber	EA.	9
MC-CG001	C-Series AIO Baseplate Glass Fiber MR	EA.	9
RC-EBF01	Extension Board Carbon Fiber	EA.	9
RM-EBG01	Extension Board Glass Fiber MR	EA.	9
MC-CF002	C-Series Head-Shoulder-Breast Baseplate Carbon Fiber	EA.	10
MC-CG002	C-Series Head-Shoulder-Breast Baseplate Glass Fiber MR	EA.	10
MC-CF004	C-Series Head-Shoulder Baseplate Carbon Fiber	EA.	10
MC-CG004	C-Series Head-Shoulder Baseplate Glass Fiber MR	EA.	10
MC-CF003	C-Series Pelvic Baseplate Carbon Fiber	EA.	10
MC-CG003	C-Series Pelvic Baseplate Glass Fiber MR	EA.	10
MC-CG001-A01	Pin-Lock Shoulder Retractor	EA.	11
MC-CG002-A02	Pin-Lock Shoulder Retractor	EA.	11
MC-CG001-A04	Leg Separator	EA.	11
MC-CG001-A03	Groin Lock	EA.	11
RC-ASF01	Overhead Arm Support Carbon Fiber	EA.	13
RM-ASG01	Overhead Arm Support Glass Fiber MR	EA.	13

Supine Breast Cushions

REF	Product Description	Pkg	Page
MC-RXE001	Supine Breast Cushions	SET	18
MC-RXE001-A01	Wedge Cushion	EA.	18
MC-RXE001-A02	Wedge Cushion	EA.	18
MC-RXE001-A03	Arm Support Cushion	EA.	18
MC-RXE001-A04	Arm Support Cushion	EA.	18

Meicen A-Series Baseplates

REF	Product Description	Pkg	Page
MC-AF001	A-Series AIO Baseplate Carbon Fiber	EA.	20
MC-AG001	A-Series AIO Baseplate Glass Fiber MR	EA.	20
MC-AG001-A01	Shoulder Retractor	EA.	20
MC-AF002	A-Series Head-Shoulder-Breast Baseplate Carbon Fiber	EA.	21
MC-AG002	A-Series Head-Shoulder-Breast Baseplate Glass Fiber MR	EA.	21
MC-AF003	A-Series Pelvic Baseplate Carbon Fiber	EA.	21
MC-AG003	A-Series Pelvic Baseplate Glass Fiber MR	EA.	21
MC-AF004	A-Series Head-Shoulder Baseplate Carbon Fiber	EA.	21
MC-AG004	A-Series Head-Shoulder Baseplate Glass Fiber MR	EA.	21

Purchasing Index

Meicen B-Series Baseplates

REF	Product Description	Pkg	Page
MC-BF003	B-Series U-Frame Head Baseplate Carbon Fiber	EA.	23
MC-BG003	B-Series U-Frame Head Baseplate Glass Fiber MR	EA.	23
MC-BF004	B-Series Tilting U-Frame Head Baseplate Carbon Fiber	EA.	23
MC-BG004	B-Series Tilting U-Frame Head Baseplate Glass Fiber MR	EA.	23

Meicen Wing Board & Supine Breast Board

REF	Product Description	Pkg	Page
MC-WBF001	Wing Board T Grip Carbon Fiber	EA.	24
MC-WBF002	Wing Board U Grip Carbon Fiber	EA.	24
MC-WBG001	Wing Board T Grip Glass Fiber MR	EA.	24
MC-WBF001-A01	T-Grip, Carbon Fiber	EA.	24
MC-WBG001-A01	T-Grip, MR	EA.	24
MC-RXF002	Supine Breast Board Carbon Fiber	EA.	25
MC-RXG002	Supine Breast Board Glass Fiber MR	EA.	25

Meicen Prone Breast & Belly Solutions

REF	Product Description	Pkg	Page
MC-RBF001	Prone Breast Board Carbon Fiber	SET	26
MC-RBG001	Prone Breast Board Glass Fiber MR	SET	26
M2422-VHT289W	Trunk Mask, 22% perforation, 2.4mm	10/BOX	26
M3038-VHT550W	Trunk Mask, 38% perforation, 3.0mm	10/BOX	26
MC-RBE001	Prone Breast Cushion	SET	27
MC-RPE002	Prone Belly Cushion	SET	27

U-frame Violet Thermoplastic Masks for Slim Frame

REF	Product Description	Pkg	Page
M2422-VUA192W	U-Shaped Head Mask, 22% perforation, 2.4mm	10/BOX	30
M3038-VUA481W	U-Shaped Head Mask, 38% perforation, 3.0mm	10/BOX	30
M2415-VUA150W	U-Shaped Head Mask, 15% perforation, 2.4mm	10/BOX	30
M24RA-VUA580W	U-Shaped Head Mask, Reinforced, 2.4mm	10/BOX	30
M30RA-VUA610W	U-Shaped Head Mask, Reinforced, 3.0mm	10/BOX	30
M24RB-VUA640W	U-Shaped Head Mask, Reinforced, Openface, 2.4mm	10/BOX	30
M30RB-VUA670W	U-Shaped Head Mask, Reinforced, Openface, 3.0mm	10/BOX	30
M24RB-VUA641W	U-Shaped Head Mask, Reinforced, Open eyes & mouth, 2.4mm	10/BOX	30
M30RB-VUA671W	U-Shaped Head Mask, Reinforced, Open eyes & mouth, 3.0mm	10/BOX	30
M2422-VUA197W	U-Shaped Head Mask, Extended 5cm, 22% perforation, 2.4mm	10/BOX	30
M3038-VUA482W	U-Shaped Head Mask, Extended 5cm, 38% perforation, 3.0mm	10/BOX	30
M2422-VUA200W	U-Shaped Head Mask, Extended 9cm, 22% perforation, 2.4mm	10/BOX	30
M3038-VUA483W	U-Shaped Head Mask, Extended 9cm, 38% perforation, 3.0mm	10/BOX	30

U-frame Violet Thermoplastic Masks for Standard Frame

REF	Product Description	Pkg	Page
M2422-VUA190W	U-Shaped Head Mask, 22% perforation, 2.4mm	10/BOX	31
M3038-VUA480W	U-Shaped Head Mask, 38% perforation, 3.0mm	10/BOX	31
M24RA-VUA581W	U-Shaped Head Mask, Reinforced, 2.4mm	10/BOX	31
M30RA-VUA611W	U-Shaped Head Mask, Reinforced, 3.0mm	10/BOX	31
M24RB-VUA642W	U-Shaped Head Mask, Reinforced, Openface, 2.4mm	10/BOX	31
M30RB-VUA672W	U-Shaped Head Mask, Reinforced, Openface, 3.0mm	10/BOX	31
M24RB-VUA643W	U-Shaped Head Mask, Reinforced, Open eyes & mouth, 2.4mm	10/BOX	31
M30RB-VUA673W	U-Shaped Head Mask, Reinforced, Open eyes & mouth, 3.0mm	10/BOX	31
M2422-VUA194W	U-Shaped Head Mask, Extended 5cm, 22% perforation, 2.4mm	10/BOX	31
M3038-VUA484W	U-Shaped Head Mask, Extended 5cm, 38% perforation, 3.0mm	10/BOX	31
M2422-VUA202W	U-Shaped Head Mask, Extended 9cm, 22% perforation, 2.4mm	10/BOX	31
M3038-VUA485W	U-Shaped Head Mask, Extended 9cm, 38% perforation, 3.0mm	10/BOX	31

Purchasing Index

S-type Violet Thermoplastic Masks

REF	Product Description	Pkg	Page
M2422-VSA190W	S-Shaped Head Mask, 22% perforation, 2.4mm	10/BOX	32
M3038-VSA480W	S-Shaped Head Mask, 38% perforation, 3.0mm	10/BOX	32
M2415-VSA150W	S-Shaped Head Mask, 15% perforation, 2.4mm	10/BOX	32
M24RA-VSA580W	S-Shaped Head Mask, Reinforced, 2.4mm	10/BOX	32
M30RA-VSA610W	S-Shaped Head Mask, Reinforced, 3.0mm	10/BOX	32
M2422-VSA192W	S-Shaped Head Mask, Extended 5cm, 22% perforation, 2.4mm	10/BOX	32
M3038-VSA482W	S-Shaped Head Mask, Extended 5cm, 38% perforation, 3.0mm	10/BOX	32
M2422-VSA193W	S-Shaped Head Mask, Extended 9cm, 22% perforation, 2.4mm	10/BOX	32
M3038-VSA483W	S-Shaped Head Mask, Extended 9cm, 38% perforation, 3.0mm	10/BOX	32
M24RA-VSA581W	S-Shaped Head Mask, Reinforced, Extended 5cm, 2.4mm	10/BOX	32
M30RA-VSA611W	S-Shaped Head Mask, Reinforced, Extended 5cm, 3.0mm	10/BOX	32
M24RA-VSA582W	S-Shaped Head Mask, Reinforced, Extended 9cm, 2.4mm	10/BOX	32
M30RA-VSA612W	S-Shaped Head Mask, Reinforced, Extended 9cm, 3.0mm	10/BOX	32
M2422-VSB190W	S-Shaped Head-Neck-Shoulder Mask, 22% perforation, 2.4mm	10/BOX	33
M3038-VSB480W	S-Shaped Head-Neck-Shoulder Mask, 38% perforation, 3.0mm	10/BOX	33
M2415-VSB150W	S-Shaped Head-Neck-Shoulder Mask, 15% perforation, 2.4mm	10/BOX	33
M24RA-VSB580W	S-Shaped Head-Neck-Shoulder Mask, Reinforced, 2.4mm	10/BOX	33
M30RA-VSB610W	S-Shaped Head-Neck-Shoulder Mask, Reinforced, 3.0mm	10/BOX	33
M2422-VSB191W	S-Shaped Head-Neck-Shoulder Mask, 22% perforation, 2.4mm, narrow shoulders	10/BOX	33
M3038-VSB481W	S-Shaped Head-Neck-Shoulder Mask, 38% perforation, 3.0mm, narrow shoulders	10/BOX	33
M24RA-VSB581W	S-Shaped Head-Neck-Shoulder Mask, Reinforced, 2.4mm, narrow shoulders	10/BOX	33
M30RA-VSB611W	S-Shaped Head-Neck-Shoulder Mask, Reinforced, 3.0mm, narrow shoulders	10/BOX	33
M24RB-VSA643W	S-Shaped Head Mask, Reinforced, Open eyes & mouth, 2.4mm	10/BOX	34
M30RB-VSA672W	S-Shaped Head Mask, Reinforced, Open eyes & mouth, 3.0mm	10/BOX	34
M24RB-VSA640W	S-Shaped Head Mask, Reinforced, Openface, 2.4mm	10/BOX	34
M30RB-VSA670W	S-Shaped Head Mask, Reinforced, Openface, 3.0mm	10/BOX	34
M24RB-VSA641W	S-Shaped Head Mask, Reinforced, Openface, 2.4mm	10/BOX	34
M30RB-VSA671W	S-Shaped Head Mask, Reinforced, Openface, 3.0mm	10/BOX	34
M24RB-VSB642W	S-Shaped Head-Neck-Shoulder Mask, Reinforced, Openface, 2.4mm	10/BOX	35
M30RB-VSB670W	S-Shaped Head-Neck-Shoulder Mask, Reinforced, Openface, 3.0mm	10/BOX	35
M24RB-VSB643W	S-Shaped Head-Neck-Shoulder Mask, Reinforced, Open eyes & mouth, 2.4mm	10/BOX	35
M30RB-VSB672W	S-Shaped Head-Neck-Shoulder Mask, Reinforced, Open eyes & mouth, 3.0mm	10/BOX	35
M24RB-VSB641W	S-Shaped Head-Neck-Shoulder Mask, Reinforced, Openface, 2.4mm	10/BOX	35
M30RB-VSB671W	S-Shaped Head-Neck-Shoulder Mask, Reinforced, Openface, 3.0mm	10/BOX	35
M24RB-VSB644W	S-Shaped Head-Neck-Shoulder Mask, Reinforced, Open face & neck, 2.4mm	10/BOX	35
M30RB-VSB673W	S-Shaped Head-Neck-Shoulder Mask, Reinforced, Open face & neck, 3.0mm	10/BOX	35

L-type Violet Thermoplastic Masks

REF	Product Description	Pkg	Page
M2422-VLA191W	Triangular Head Mask, 22% perforation, with nose hole, 2.4mm	10/BOX	36
M3038-VLA480W	Triangular Head Mask, 38% perforation, with nose hole, 3.0mm	10/BOX	36
M2438-VLA320W	Triangular Head Mask, 38% perforation, with nose hole, 2.4mm, for pediatric	10/BOX	36
M2422-VLA192W	Triangular Head Mask, Extended 3cm, 22% perforation, with nose hole, 2.4mm	10/BOX	36
M3038-VLA481W	Triangular Head Mask, Extended 3cm, 38% perforation, with nose hole, 3.0mm	10/BOX	36
M2415-VLA151W	Triangular Head Mask, Extended 3cm, 15% perforation, with nose hole, 2.4mm	10/BOX	36
M24RA-VLA580W	Triangular Head Mask, Reinforced, with nose hole, 2.4mm	10/BOX	36
M30RA-VLA610W	Triangular Head Mask, Reinforced, with nose hole, 3.0mm	10/BOX	36
M2422-VLE190W	Neck-Shoulder Mask, 22% perforation, with nose hole, 2.4mm	10/BOX	36
M3038-VLE480W	Neck-Shoulder Mask, 38% perforation, with nose hole, 3.0mm	10/BOX	36
M2422-VLE191W	Neck-Shoulder Mask, 22% perforation, with nose hole, 2.4mm, for pediatric	10/BOX	36
M2422-VLF192W	Neck-Thorax Mask, 22% perforation, 2.4mm	10/BOX	36
M3038-VLF480W	Neck-Thorax Mask, 38% perforation, 3.0mm	10/BOX	36
M2422-VLB191W	Head-Neck-Shoulder Mask, 22% perforation, with nose hole, 2.4mm	10/BOX	37
M3038-VLB481W	Head-Neck-Shoulder Mask, 38% perforation, with nose hole, 3.0mm	10/BOX	37
M2422-VLB192W	Head-Neck-Shoulder Mask, 22% perforation, with nose hole, 2.4mm, for pediatric	10/BOX	37
M24RA-VLB580W	Head-Neck-Shoulder Mask, Reinforced, with nose hole, 2.4mm	10/BOX	37
M30RA-VLB610W	Head-Neck-Shoulder Mask, Reinforced, with nose hole, 3.0mm	10/BOX	37
M3038-VLK483W	Chest-Pelvic Mask, 38% perforation, 3.0mm, 228 clamp	10/BOX	37
M3038-VLK482W	Chest-Pelvic Mask, 38% perforation, 3.0mm, 243 clamp	10/BOX	37
M3038-VLK481W	Chest-Pelvic Mask, 38% perforation, 3.0mm, 247 clamp	10/BOX	37
M3038-VLL483W	Pelvic Mask, 38% perforation, 3.0mm, 228 clamp	10/BOX	37
M3038-VLL482W	Pelvic Mask, 38% perforation, 3.0mm, 243 clamp	10/BOX	37
M3038-VLL481W	Pelvic Mask, 38% perforation, 3.0mm, 247 clamp	10/BOX	37
M2422-VLG190W	Thorax-Abdomen Mask, 22% perforation, 2.4mm	10/BOX	37
M3038-VLG480W	Thorax-Abdomen Mask, 38% perforation, 3.0mm	10/BOX	37
M24RB-VLA641W	Triangular Head Mask, Reinforced, Openface, 2.4mm	10/BOX	38
M30RB-VLA671W	Triangular Head Mask, Reinforced, Openface, 3.0mm	10/BOX	38
M24RB-VLA643W	Triangular Head Mask, Reinforced, Open eyes & mouth, 2.4mm	10/BOX	38
M30RB-VLA673W	Triangular Head Mask, Reinforced, Open eyes & mouth, 3.0mm	10/BOX	38
M24RB-VLB642W	Head-Neck-Shoulder Mask, Reinforced, Open eyes & mouth, 2.4mm	10/BOX	38
M30RB-VLB672W	Head-Neck-Shoulder Mask, Reinforced, Open eyes & mouth, 3.0mm	10/BOX	38
M24RB-VLB643W	Head-Neck-Shoulder Mask, Openface, 2.4mm	10/BOX	38
M30RB-VLB673W	Head-Neck-Shoulder Mask, Openface, 3.0mm	10/BOX	38
M24RB-VLC641W	Head-Neck-Shoulder Mask, Reinforced, Extended, Openface, 2.4mm	10/BOX	38
M30RB-VLC671W	Head-Neck-Shoulder Mask, Reinforced, Extended, Openface, 3.0mm	10/BOX	38

Purchasing Index

P-type Violet Thermoplastic Masks

REF	Product Description	Pkg	Page
M2422-VPA190W	Triangular Head Mask, 22% perforation, with nose hole, 2.4mm	10/BOX	39
M3038-VPA480W	Triangular Head Mask, 38% perforation, with nose hole, 3.0mm	10/BOX	39
M2415-VPA150W	Triangular Head Mask, 15% perforation, with nose hole, 2.4mm	10/BOX	39
M24RA-VPA580W	Triangular Head Mask, Reinforced, with nose hole, 2.4mm	10/BOX	39
M30RA-VPA610W	Triangular Head Mask, Reinforced, with nose hole, 3.0mm	10/BOX	39
M2422-VPA191W	Triangular Head Mask, Extended 5cm, 22% perforation, with nose hole, 2.4mm	10/BOX	39
M3038-VPA481W	Triangular Head Mask, Extended 5cm, 38% perforation, with nose hole, 3.0mm	10/BOX	39
M2415-VPA151W	Triangular Head Mask, Extended 5cm, 15% perforation, with nose hole, 2.4mm	10/BOX	39
M2422-VPE190W	Neck-Shoulder Mask, 22% perforation, with nose hole, 2.4mm	10/BOX	39
M3038-VPE480W	Neck-Shoulder Mask, 38% perforation, with nose hole, 3.0mm	10/BOX	39
M2422-VPB190W	Head-Neck-Shoulder Mask, 22% perforation, with nose hole, 2.4mm	10/BOX	39
M3038-VPB480W	Head-Neck-Shoulder Mask, 38% perforation, with nose hole, 3.0mm	10/BOX	39
M2415-VPB150W	Head-Neck-Shoulder Mask, 15% perforation, with nose hole, 2.4mm	10/BOX	39
M24RA-VPB580W	Head-Neck-Shoulder Mask, Reinforced, with nose hole, 2.4mm	10/BOX	40
M30RA-VPB610W	Head-Neck-Shoulder Mask, Reinforced, with nose hole, 3.0mm	10/BOX	40
M2422-VPK190W	Chest-Pelvic Mask, 22% perforation, 2.4mm	10/BOX	40
M3038-VPK481W	Chest-Pelvic Mask, 38% perforation, 3.0mm	10/BOX	40
M2422-VPA192W	Triangular Head Mask, 22% perforation, with nose hole, 2.4mm, for pediatric	10/BOX	40
M3038-VPA482W	Triangular Head Mask, 38% perforation, with nose hole, 3.0mm, for pediatric	10/BOX	40
M2422-VPB191W	Head-Neck-Shoulder Mask, 22% perforation, with nose hole, 2.4mm, for pediatric	10/BOX	40
M3038-VPB481W	Head-Neck-Shoulder Mask, 38% perforation, with nose hole, 3.0mm, for pediatric	10/BOX	40
M24RA-VPB581W	Head-Neck-Shoulder Mask, Reinforced, with nose hole, 2.4mm, for pediatric	10/BOX	40
M30RA-VPB611W	Head-Neck-Shoulder Mask, Reinforced, with nose hole, 3.0mm, for pediatric	10/BOX	40
M24RB-VPA640W	Triangular Head Mask, Openface, 2.4mm	10/BOX	41
M30RB-VPA670W	Triangular Head Mask, Openface, 3.0mm	10/BOX	41
M24RB-VPA641W	Triangular Head Mask, Reinforced, Openface, 2.4mm	10/BOX	41
M30RB-VPA671W	Triangular Head Mask, Reinforced, Openface, 3.0mm	10/BOX	41
M24RB-VPA642W	Triangular Head Mask, Reinforced, Open eyes & mouth, 2.4mm	10/BOX	41
M30RB-VPA672W	Triangular Head Mask, Reinforced, Open eyes & mouth, 3.0mm	10/BOX	41
M24RB-VPB640W	Head-Neck-Shoulder Mask, Reinforced, Openface, 2.4mm	10/BOX	41
M30RB-VPB670W	Head-Neck-Shoulder Mask, Reinforced, Openface, 3.0mm	10/BOX	41
M24RB-VPB641W	Head-Neck-Shoulder Mask, Reinforced, Open eyes & mouth, 2.4mm	10/BOX	41
M30RB-VPB671W	Head-Neck-Shoulder Mask, Reinforced, Open eyes & mouth, 3.0mm	10/BOX	41

Pin-lock type Violet Thermoplastic Masks

REF	Product Description	Pkg	Page
M24RA-VSF580W	Neck-Thorax Mask, Reinforced, 2.4mm	10/BOX	46
M30RA-VSF610W	Neck-Thorax Mask, Reinforced, 3.0mm	10/BOX	46
M2422-VST190W	Trunk Mask (Chest-Pelvic), Pin-Lock, 22% perforation, 2.4mm	10/BOX	46
M3038-VST480W	Trunk Mask (Chest-Pelvic), Pin-Lock, 38% perforation, 3.0mm	10/BOX	46
M2422-VST191W	Trunk Mask (Pelvic), Pin-Lock, with Groin Lock, 22% perforation, 2.4mm	10/BOX	46
M3038-VST481W	Trunk Mask (Pelvic), Pin-Lock, with Groin Lock, 38% perforation, 3.0mm	10/BOX	46
M2422-VSH190W	Breast Mask, Pin-Lock, 22% perforation, 2.4mm	10/BOX	46
M3038-VSH480W	Breast Mask, Pin-Lock, 38% perforation, 3.0mm	10/BOX	46

H-type Violet Thermoplastic Masks

REF	Product Description	Pkg	Page
M2422-VHT190W	Trunk Mask (Chest-Pelvic), 22% perforation, 2.4mm	10/BOX	43
M3038-VHT480W	Trunk Mask (Chest-Pelvic), 38% perforation, 3.0mm	10/BOX	43
M2422-VHT191W	Trunk Mask (Pelvic), with Groin Lock, 22% perforation, 2.4mm	10/BOX	43
M3038-VHT481W	Trunk Mask (Pelvic), with Groin Lock, 38% perforation, 3.0mm	10/BOX	43
M24RA-VHT580W	Trunk Mask (Chest-Pelvic), Reinforced, 2.4mm	10/BOX	43
M30RA-VHT610W	Trunk Mask (Chest-Pelvic), Reinforced, 3.0mm	10/BOX	43
M2400-VHT120W	Trunk Mask (Chest-Pelvic), Plain, 2.4mm	10/BOX	43
M3000-VHT420W	Trunk Mask (Chest-Pelvic), Plain, 3.0mm	10/BOX	43

M-type Violet Thermoplastic Masks

REF	Product Description	Pkg	Page
M2422-VMA190W	Triangular Head Mask, 22% perforation, with nose hole, 2.4mm	10/BOX	44
M3038-VMA480W	Triangular Head Mask, 38% perforation, with nose hole, 3.0mm	10/BOX	44
M24RA-VMA580W	Triangular Head Mask, Reinforced, with nose hole, 2.4mm	10/BOX	44
M30RA-VMA610W	Triangular Head Mask, Reinforced, with nose hole, 3.0mm	10/BOX	44
M2422-VMB190W	Head-Neck-Shoulder Mask, 22% perforation, with nose hole, 2.4mm	10/BOX	44
M3038-VMB480W	Head-Neck-Shoulder Mask, 38% perforation, with nose hole, 3.0mm	10/BOX	44
M24RA-VMB580W	Head-Neck-Shoulder Mask, Reinforced, with nose hole, 2.4mm	10/BOX	44
M30RA-VMB610W	Head-Neck-Shoulder Mask, Reinforced, with nose hole, 3.0mm	10/BOX	44
M2422-VME190W	Neck-Shoulder Mask, 22% perforation, with nose hole, 2.4mm	10/BOX	44
M3038-VME480W	Neck-Shoulder Mask, 38% perforation, with nose hole, 3.0mm	10/BOX	44

Purchasing Index

X-Knife Violet Thermoplastic Masks

REF	Product Description	Pkg	Page
M3038-VXA481W	X-Knife Mask, 38% perforation, 3.0mm	5/BOX	45
M3038-VXB480W	X-Knife Head-Shoulder Mask, 38% perforation, 3.0mm	5/BOX	45
M3038-VXA490W	X-Knife Mask, 38% perforation, 3.0mm	10/BOX	45

E-type Violet Thermoplastic Masks

REF	Product Description	Pkg	Page
M24RA-VEA600W	Triangular Head Mask, Reinforced, 2.4mm	10/BOX	46
M30RA-VEA620W	Triangular Head Mask, Reinforced, 3.0mm	10/BOX	46
M24RA-VEA602W	Triangular Head Mask, Reinforced, 2.4mm	10/BOX	46
M30RA-VEA621W	Triangular Head Mask, Reinforced, 3.0mm	10/BOX	46
M2422-VEA190W	Head Mask, 22% perforation, 2.4mm	10/BOX	46
M3038-VEA480W	Head Mask, 38% perforation, 3.0mm	10/BOX	46
M2415-VEA150W	Head Mask, 15% perforation, 2.4mm	10/BOX	46
M24RA-VEA580W	Head Mask, Reinforced, 2.4mm	10/BOX	46
M30RA-VEA610W	Head Mask, Reinforced, 3.0mm	10/BOX	46
M2422-VEA192W	Head Mask, Extended 5cm, 22% perforation, 2.4mm	10/BOX	47
M3038-VEA482W	Head Mask, Extended 5cm, 38% perforation, 3.0mm	10/BOX	47
M2422-VEA193W	Head Mask, Extended 21cm, 22% perforation, 2.4mm	10/BOX	47
M3038-VEA483W	Head Mask, Extended 21cm, 38% perforation, 3.0mm	10/BOX	47
M2422-VEB190W	Head-Neck-Shoulder Mask, 22% perforation, 2.4mm	10/BOX	47
M3038-VEB480W	Head-Neck-Shoulder Mask, 38% perforation, 3.0mm	10/BOX	47
M2415-VEB150W	Head-Neck-Shoulder Mask, 15% perforation, 2.4mm	10/BOX	47
M24RA-VEB580W	Head-Neck-Shoulder Mask, Reinforced, 2.4mm	10/BOX	47
M30RA-VEB610W	Head-Neck-Shoulder Mask, Reinforced, 3.0mm	10/BOX	47
M3038-VEK480W	Chest-Pelvic Mask, 38% perforation, 3.0mm	10/BOX	47

V-type Violet Thermoplastic Masks

REF	Product Description	Pkg	Page
M2422-VVA190W	V-Shaped Head Mask, 22% perforation, 2.4mm	10/BOX	52
M3038-VVA480W	V-Shaped Head Mask, 38% perforation, 3.0mm	10/BOX	52

REF	Product Description	Pkg	Page
M24RA-VVA580W	V-Shaped Head Mask, Reinforced, 2.4mm	10/BOX	48
M30RA-VVA610W	V-Shaped Head Mask, Reinforced, 3.0mm	10/BOX	48
M24RB-VVA640W	V-Shaped Head Mask, Reinforced, Open eyes & mouth, 2.4mm	10/BOX	48
M30RB-VVA670W	V-Shaped Head Mask, Reinforced, Open eyes & mouth, 3.0mm	10/BOX	48
M24RB-VVA641W	V-Shaped Head Mask, Reinforced, Openface, 2.4mm	10/BOX	48
M30RB-VVA671W	V-Shaped Head Mask, Reinforced, Openface, 3.0mm	10/BOX	48
M2422-VVA191W	V-Shaped Head Mask, Extended 5cm, 22% perforation, 2.4mm	10/BOX	48
M3038-VVA481W	V-Shaped Head Mask, Extended 5cm, 38% perforation, 3.0mm	10/BOX	48
M2422-VVA192W	V-Shaped Head Mask, Extended 9cm, 22% perforation, 2.4mm	10/BOX	48
M3038-VVA482W	V-Shaped Head Mask, Extended 9cm, 38% perforation, 3.0mm	10/BOX	48
M2422-VVB190W	V-Shaped Head-Neck-Shoulder Mask, 22% perforation, 2.4mm	10/BOX	49
M3038-VVB480W	V-Shaped Head-Neck-Shoulder Mask, 38% perforation, 3.0mm	10/BOX	49
M24RA-VVB580W	V-Shaped Head-Neck-Shoulder Mask, Reinforced, 2.4mm	10/BOX	49
M30RA-VVB610W	V-Shaped Head-Neck-Shoulder Mask, Reinforced, 3.0mm	10/BOX	49
M24RB-VVB640W	V-Shaped Head-Neck-Shoulder Mask, Reinforced, Openface, 2.4mm	10/BOX	49
M30RB-VVB670W	V-Shaped Head-Neck-Shoulder Mask, Reinforced, Openface, 3.0mm	10/BOX	49
M24RB-VVB641W	V-Shaped Head-Neck-Shoulder Mask, Reinforced, Open eyes & mouth, 2.4mm	10/BOX	49
M30RB-VVB671W	V-Shaped Head-Neck-Shoulder Mask, Reinforced, Open eyes & mouth, 3.0mm	10/BOX	49

Meicen Vacuum Cushions

REF	Product Description	Pkg	Page
RC-E5003	100cm x 70cm, fill 30L, For Chest-Pelvic Support	EA.	51
RM-E5003	100cm x 70cm, fill 30L, For Chest-Pelvic Support, MR	EA.	51
RC-E5005S	100cm x 70cm, fill 30L, W/Indexed, For Chest-Pelvic Support	EA.	51
RM-E5005S	100cm x 70cm, fill 30L, W/Indexed, For Chest-Pelvic Support, MR	EA.	51
RC-E4013	70cm x 70cm, fill 20L, For Breast Support	EA.	51
RM-E4013	70cm x 70cm, fill 20L, For Breast Support, MR	EA.	51
RC-E4014F	70cm x 70cm, fill 20L, W/Indexed, For Breast Support	EA.	51
RM-E4014F	70cm x 70cm, fill 20L, W/Indexed, For Breast Support, MR	EA.	51
RC-E5024	120cm x 80cm, fill 32L, For Chest-Pelvic Support	EA.	51
RM-E5024	120cm x 80cm, fill 32L, For Chest-Pelvic Support, MR	EA.	51
RC-E5025S	120cm x 80cm, fill 32L, W/Indexed, For Chest-Pelvic Support	EA.	51
RM-E5025S	120cm x 80cm, fill 32L, W/Indexed, For Chest-Pelvic Support, MR	EA.	51
RC-E5010	100cm x 75cm, fill 30L, For Chest-Pelvic Support	EA.	51
RM-E5010	100cm x 75cm, fill 30L, For Chest-Pelvic Support, MR	EA.	51

Purchasing Index

REF	Product Description	Pkg	Page
RC-E5012S	100cm x 75cm, fill 30L, W/Indexed, For Chest-Pelvic Support	EA.	51
RM-E5012S	100cm x 75cm, fill 30L, W/Indexed, For Chest-Pelvic Support, MR	EA.	51
RC-E5023S	120cm x 75cm, fill 34L, W/Indexed, For Chest-Pelvic Support	EA.	51
RM-E5023S	120cm x 75cm, fill 34L, W/Indexed, For Chest-Pelvic Support, MR	EA.	51
RC-E5032	150cm x 75cm, fill 45L, For Chest-Pelvic Support	EA.	51
RM-E5032	150cm x 75cm, fill 45L, For Chest-Pelvic Support, MR	EA.	51
RC-E5033S	150cm x 75cm, fill 45L, W/Indexed, For Chest-Pelvic Support	EA.	51
RM-E5033S	150cm x 75cm, fill 45L, W/Indexed, For Chest-Pelvic Support, MR	EA.	51
RC-E1005	75cm x 75cm, fill 15L, Arms up, For Breast Support	EA.	52
RM-E1005	75cm x 75cm, fill 15L, Arms up, For Breast Support, MR	EA.	52
RC-E3005F	125cm x 110cm, fill 31L, T-Shaped, For Thorax-Abdomen Support	EA.	52
RM-E3005F	125cm x 110cm, fill 31L, T-Shaped, For Thorax-Abdomen Support, MR	EA.	52
RC-E3010F	125cm x 110cm, fill 32L, T-Shaped, For Thorax-Abdomen Support	EA.	52
RM-E3010F	125cm x 110cm, fill 32L, T-Shaped, For Thorax-Abdomen Support, MR	EA.	52
RC-E3019	105cm x 75cm, fill 28L, Arms up, For Thorax-Abdomen Support	EA.	52
RM-E3019	105cm x 75cm, fill 28L, Arms up, For Thorax-Abdomen Support, MR	EA.	52
RC-E4003	110cm x 55cm, fill 15L, Wing board breast, For Breast Support	EA.	52
RM-E4003	110cm x 55cm, fill 15L, Wing board breast, For Breast Support, MR	EA.	52
RC-E4004	110cm x 55cm, fill 19L, Wing board breast, For Breast Support	EA.	52
RM-E4004	110cm x 55cm, fill 19L, Wing board breast, For Breast Support, MR	EA.	52
RC-E1004	61cm x 65cm, fill 8L, S-Shaped, For Head-Shoulder Support	EA.	52
RM-E1004	61cm x 65cm, fill 8L, S-Shaped, For Head-Shoulder Support, MR	EA.	52

Diagnostic Couchtop Overlays

REF	Product Description	Pkg	Page
MC-CTF019	CT Overlay	EA.	57
MC-MRG002	MRI Overlay	EA.	57

Accessories

REF	Product Description	Pkg	Page
RE-HSF6	Cabon Fiber Headrest, set of 6	EA.	54
RE-HSB5	Color Headrest, S-type, set of 5	EA.	54
RE-HSC6	Clear Headrest, set of 6	EA.	54
RE-HSP5	Color Headrest, A-type, set of 5	EA.	54
RE-HSA04	Shim for Head	EA.	54
RE-HSA05	Shim for Head-Neck-Shoulder	EA.	54
RC-LBE01	For Elekta Couchtops, Aluminum	EA.	55
RC-LBE02	For Elekta Couchtops, Carbon Fiber	EA.	55
RM-LBE01	For Elekta Couchtops, Glass Fiber MR	EA.	55
RC-LBV01	For Varian Couchtops, Aluminum	EA.	55
RC-LBV02	For Varian Couchtops, Carbon Fiber	EA.	55
RM-LBV01	For Varian Couchtops, Glass Fiber MR	EA.	55
RC-LBC07	For Vacuum Cushions, Carbon Fiber	EA.	55
RC-LBV03	For Meicen C-Series AIO Baseplate & Varian Couchtops, Carbon Fiber	EA.	55
RC-LBE03	For Meicen C-Series AIO Baseplate & Elekta Couchtops, Carbon Fiber	EA.	55
RC-LBC02	Vacuum Cushion Positioners Strip, Carbon Fiber	EA.	55
R-PM002	Vauum Pump	EA.	55
MC-CG001-A01	Pin-Lock Shoulder Retractor	EA.	56
MC-CG002-A02	Pin-Lock Shoulder Retractor	EA.	56
MC-YG001-A03	Pin-Lock Shoulder Retractor	EA.	56
MC-AG001-A01	Shoulder Retractor	EA.	56
RC-HGF01	Hand Grips, Carbon Fiber	EA.	56
RM-HGG01	Hand Grips, MR	EA.	56
RE-FCA01	Knee Cushion	EA.	57
RE-FCA02	Elevation Cushion	EA.	57
RE-FCA03	Feet Cushion	EA.	57
04RF-G2	Water Bath	EA.	57
04126	Bolus 300x300x5mm	EA.	57
04127	Bolus 300x300x10mm	EA.	57