

Meicen S-Series SBRT Solution

Meicen S-Series SBRT Solution is highly accurate, comfortable, multi-functional and user-friendly device for SBRT treatment, which is compatible with the couchtops of all brands. The paddle design on the belly and knee bridges can be adjusted in the direction of left and right, up and down to control the breathing range of the patient effectively during the treatment. The operation of quick-lock and unlock system in bridges is effortless. And, there are scales to record the fixed positions.

The solution features a flexible and fully adjustable breathing bridge to suit different patient sizes. And the abundant scale marks improve the accuracy of positioning repeatability.

It can be used with Wing Board, Vacuum Cushions, Thermoplastic Masks, Supine Breast Cushion and other products to provide various positioning solutions for Radiation Therapy.

Specification: 2000(L) x 530(W)/mm



Meicen S-Series SBRT Solution



- High accuracy**
Clear and precise scale, improves the accuracy of positioning repeatability.
- Versatility**
Multifunctional design, simple to form a set-up for various demand.
- Patient comfort**
Ergonomically design for patient comfort.
- User-friendly**
Easy and simple to operate.




Meicen S-Series SBRT Solution Components

SBRT Platform

MC-STF001-A01 Carbon Fiber

MC-STG001-A01 Glass Fiber MR

- Multiple corresponding holes, efficient indexing, install without restraint.
- Lightweight, easy to move.

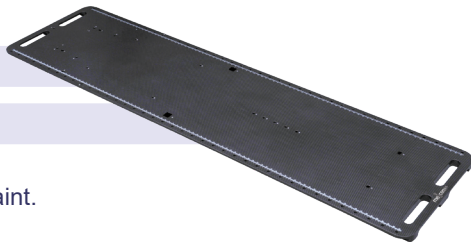

Wing Board T-Grip (Page 23)

MC-WBF001 Carbon Fiber

MC-WBG001 Glass Fiber MR

- T-Grip design provides comfortable arms extension.
- Can be used with the T-Shaped Vacuum Cushion, RC-E3005F or RC-E3010F.



Belly Bridge

MC-STF001-A02 Carbon Fiber

MC-STG001-A02 Glass Fiber MR

- Ergonomic Compression Paddle for better comfort.
- The paddle device can be adjusted in the direction of left and right, up and down, and has scales to record its fixed position.
- Abundant scale marks improve the accuracy of positioning repeatability.



Knee Bridge

MC-STF001-A03 Carbon Fiber

MC-STG001-A03 Glass Fiber MR

- The knee paddle device can be adjusted left and right, up and down, and has scales to record its fixed position.
- Abundant scale marks improve the accuracy of positioning repeatability.
- Knee paddle device is optional in both legs or single leg immobilization.




Meicen S-Series SBRT Solution Components

Respiratory Belt

MC-STF001-A04 Carbon Fiber

MC-STG001-A04 Glass Fiber MR

- Provide pneumatic compression and immobilization to the abdominal region.

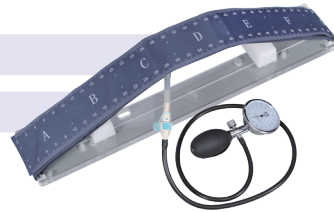

T-Shaped Vacuum Cushion (Page 58)

RC-E3005F Non-open headrest

RC-E3010F Open headrest

- Specification: 125x110cm.
- Open headrest area can be used together with the headrest to reduce the ray blocking.
- Provides arms-up positioning and support the torso region.




Knee & Feet Cushion

RE-FCA01 Knee Cushion

RE-FCA03 Feet Cushion

RE-FCA02 Elevation Cushion

- Easy to operate, improve the comfort and stability of supine position for leg positioning during fixation.
- It can be directly fixed on the couchtop for use.

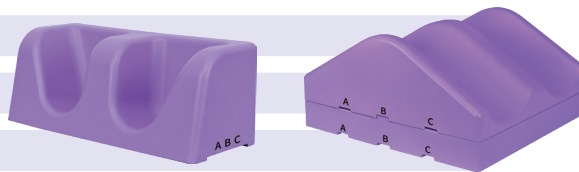

Feet Cushion


Knee & Elevation Cushion

SBRT Indexing Bar

MC-STF001-A07 Carbon Fiber

MC-STG001-A07 Glass Fiber MR

- Allows positioning devices to be fixed to the SBRT Platform for enhancing reproducibility.

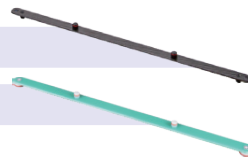

Meicen S-Series SBRT Solution



Carbon Fiber SBRT Solution Package

	Component	REF No.	Qty
MC-STF001 Standard Package	SBRT Platform	MC-STF001-A01	1
	Wing Board T-Grip	MC-WBF001	1
	Belly Bridge	MC-STF001-A02	1
	Knee Bridge	MC-STF001-A03	1
	Knee Cushion	RE-FCA01	1
	Elevation Cushion	RE-FCA02	1
	Feet Cushion	RE-FCA03	1
	T-Shaped Vacuum Cushion	RC-E3010F	1
	SBRT Indexing Bar	MC-STF001-A07	4
	Carbon Fiber Headrest	RE-HSF6-B	1
Optional Product	Carbon Fiber Headrest	RE-HSF6	1
	Color Headrest	RE-HSB5	1
	Clear Headrest	RE-HSC6	1
	Respiratory Belt	MC-STF001-A04	1
	Indexing Bar (for E-style indexed couchtops)	RC-LBE02	2
	Indexing Bar (for V-style indexed couchtops)	RC-LBV02	2
	Shoulder Retractor	MC-STG001-A08	1

Meicen S-Series SBRT Solution



MR SBRT Solution Package

	Component	REF No.	Qty
MC-STG001 Standard Package	SBRT Platform	MC-STG001-A01	1
	Wing Board T-Grip	MC-WBG001	1
	Belly Bridge	MC-STG001-A02	1
	Knee Bridge	MC-STG001-A03	1
	Knee Cushion	RE-FCA01	1
	Elevation Cushion	RE-FCA02	1
	Feet Cushion	RE-FCA03	1
	T-Shaped Vacuum Cushion	RM-E3010F	1
	SBRT Indexing Bars	MC-STG001-A07	4
	Clear Headrest	RE-HSC6-B	1
Optional Product	Color Headrest	RE-HSB5	1
	Clear Headrest	RE-HSC6	1
	Respiratory Belt	MC-STG001-A04	1
	Indexing Bar (for E-style indexed couchtops)	RM-LBE01	2
	Indexing Bar (for V-style indexed couchtops)	RM-LBV01	2
	Shoulder Retractor	MC-STG001-A08	1