


Violet Thermoplastic Masks



Meicen Violet Thermoplastic Masks

Meicen Violet Thermoplastic Masks are composed of a new, premium thermoplastic material designed specifically to offer the ultimate in patient comfort with exceptional rigidity and extremely low shrinkage, and primary focus is the patient experience. It is available in multiple perforation styles and compatible with most manufacturers' baseplates.

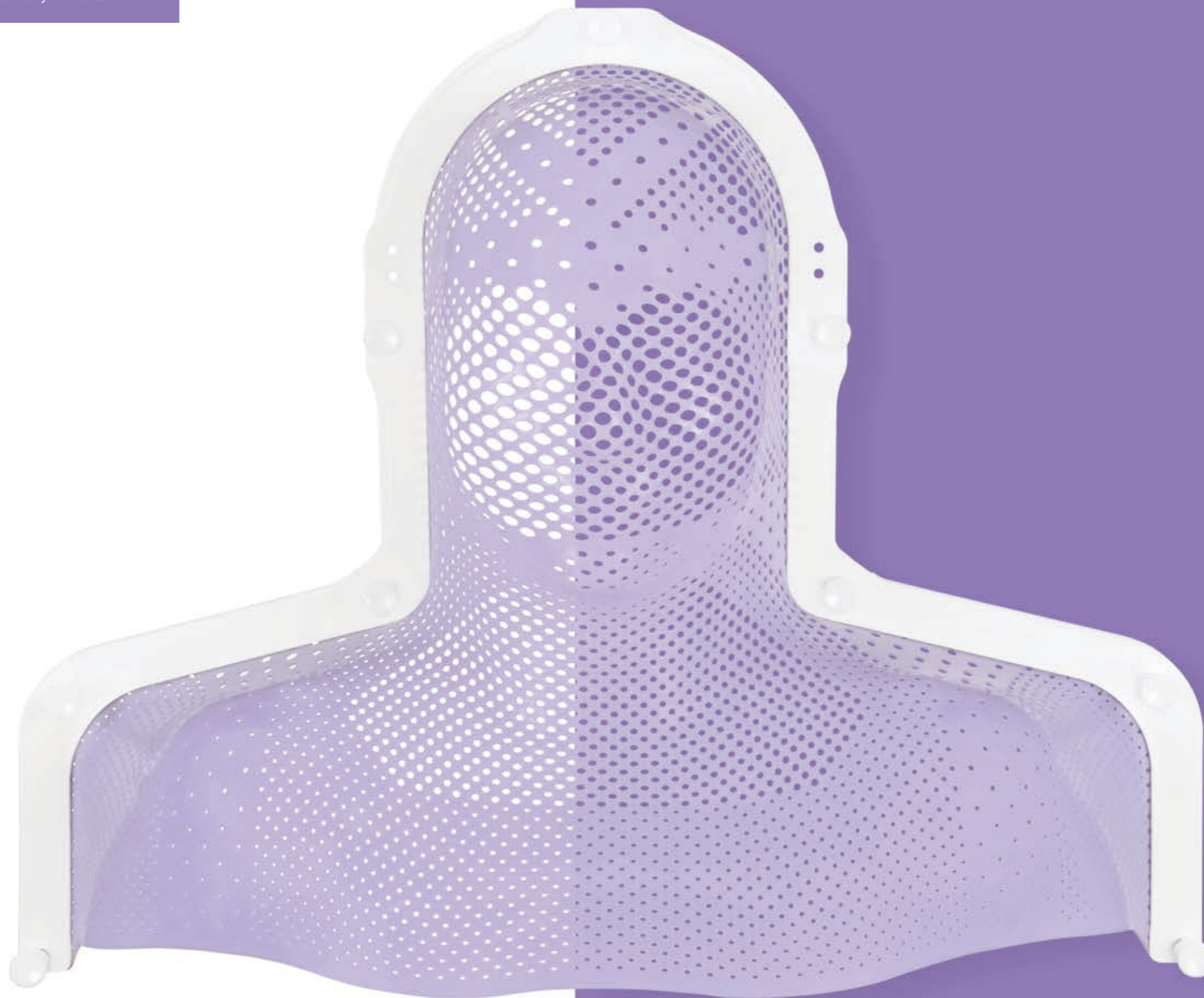

- Ergonomically designed for patient comfort
- Violet color to alleviate patient anxiety
- Optimal molding time
- Low shrinkage
- Non-Stick

Violet a traditional Chinese color that evokes tranquility and harmony. Reminiscent of subtle hues reflected in snow, it creates an ethereal, serene atmosphere. With its timeless quality, violet inspires the imagination and reflects a static beauty—calm, peaceful, and unhurried. In Chinese culture, such colors carry deep artistic significance, symbolizing the pursuit of inner peace and harmony with nature.



Meicen Violet Thermoplastic Masks


Meicen Reinforced Series Masks feature a combination of enhanced perforation and solid thermoplastic to provide rigidity and secure immobilization for IMRT treatments. Openface Masks can be compatible with SGRT (Surface Guided Radiation Therapy) camera systems.




S-type Reinforced Openface



P-type Reinforced

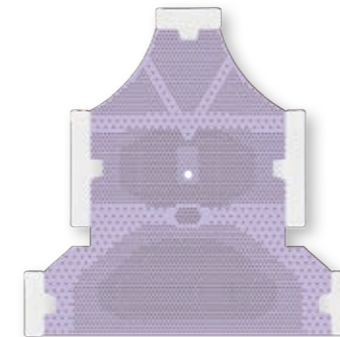

L-type Reinforced

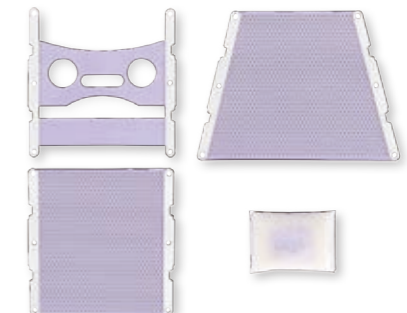

Stereotactic X-Knife

Meicen Violet Thermoplastic Masks are easy to prepare and can be heated by water bath at 65 - 70°C (149 - 158°F) or oven at 70 - 75°C (158 - 167°F).

Moulding Time: **1 - 2 minutes**

Cooling Time: **8 - 15 minutes**

Thickness 2.4mm, 3.0mm are available.

Please contact our sales representative for more information.



U-type Violet Thermoplastic Masks **Slim Frame**



U-type Violet Thermoplastic Masks **Standard Frame**



U-Shaped Head Mask

Standard

M2422-VUA192W	22% perforation, 2.4mm
M3038-VUA481W	38% perforation, 3.0mm
M2415-VUA150W	15% perforation, 2.4mm



U-Shaped Head Mask

Standard

M2422-VUA190W	22% perforation, 2.4mm
M3038-VUA480W	38% perforation, 3.0mm



U-Shaped Head Mask

Reinforced

M24RA-VUA580W	2.4mm
M30RA-VUA610W	3.0mm



U-Shaped Head Mask

Reinforced

M24RA-VUA581W	2.4mm
M30RA-VUA611W	3.0mm



U-Shaped Head Mask

Reinforced, Openface

M24RB-VUA640W	2.4mm
M30RB-VUA670W	3.0mm



U-Shaped Head Mask

Reinforced, Openface

M24RB-VUA642W	2.4mm
M30RB-VUA672W	3.0mm



U-Shaped Head Mask

Reinforced, Open eyes & mouth

M24RB-VUA641W	2.4mm
M30RB-VUA671W	3.0mm



U-Shaped Head Mask

Reinforced, Open eyes & mouth

M24RB-VUA643W	2.4mm
M30RB-VUA673W	3.0mm



U-Shaped Head Mask

Standard, Extended

Extended 5cm		Extended 9cm
M2422-VUA197W	M2422-VUA200W	22% perforation, 2.4mm
M3038-VUA482W	M3038-VUA483W	38% perforation, 3.0mm



U-Shaped Head Mask

Standard, Extended

Extended 5cm		Extended 9cm
M2422-VUA194W	M2422-VUA202W	22% perforation, 2.4mm
M3038-VUA484W	M3038-VUA485W	38% perforation, 3.0mm

Meicen U-type Masks fit all standard U-type Baseplates from Meicen and other manufacturers.

2.0mm & 3.2mm thickness are available also, please specify thickness when requesting a quote.

S-type Violet Thermoplastic Masks



S-Shaped Head Mask

Standard

M2422-VSA190W	22% perforation, 2.4mm
M3038-VSA480W	38% perforation, 3.0mm
M2415-VSA150W	15% perforation, 2.4mm



S-Shaped Head Mask

Reinforced

M24RA-VSA580W	2.4mm
M30RA-VSA610W	3.0mm



S-Shaped Head Mask

Standard, Extended

Extended 5cm	Extended 9cm	
M2422-VSA192W	M2422-VSA193W	22% perforation, 2.4mm
M3038-VSA482W	M3038-VSA483W	38% perforation, 3.0mm



S-Shaped Head Mask

Reinforced, Extended

Extended 5cm	Extended 9cm	
M24RA-VSA581W	M24RA-VSA582W	2.4mm
M30RA-VSA611W	M30RA-VSA612W	3.0mm


S-type Violet Thermoplastic Masks



S-Shaped Head-Neck-Shoulder Mask

Standard

M2422-VSB190W	22% perforation, 2.4mm
M3038-VSB480W	38% perforation, 3.0mm
M2415-VSB150W	15% perforation, 2.4mm

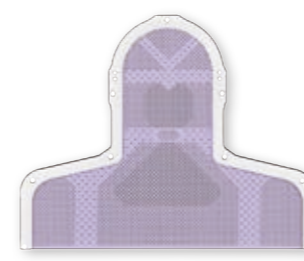


S-Shaped Head-Neck-Shoulder Mask

Reinforced

M24RA-VSB580W	2.4mm
M30RA-VSB610W	3.0mm


Specific design for Meicen C & Y-Series Baseplates

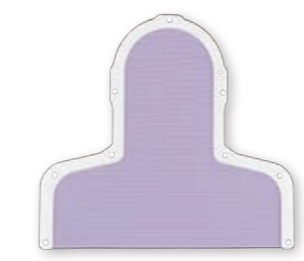

S-Shaped Head-Neck-Shoulder Mask

Standard, narrow shoulder

M2422-VSB191W	22% perforation, 2.4mm
M3038-VSB481W	38% perforation, 3.0mm

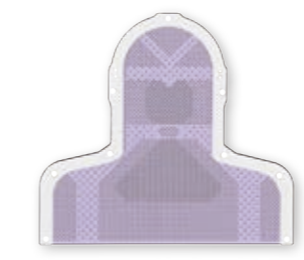

S-Shaped Head-Neck-Shoulder Mask

Reinforced, narrow shoulder

M24RA-VSB581W	2.4mm
M30RA-VSB611W	3.0mm

Meicen S-type Masks fit all standard S-type Baseplates from Meicen and other manufacturers.

2.0mm & 3.2mm thickness are available also, please specify thickness when requesting a quote.

S-type Violet Thermoplastic Masks Openface

Meicen Openface Series Masks aims to alleviate symptoms of claustrophobia during radiotherapy. All types of masks are available as openface. Adjustable pulling of the open part when apply to different patients assures comfort and reinforced effect.




S-Shaped Head Mask
Reinforced, Open eyes & mouth

M24RB-VSA643W	2.4mm
M30RB-VSA672W	3.0mm



S-Shaped Head Mask
Reinforced, Openface

M24RB-VSA640W	2.4mm
M30RB-VSA670W	3.0mm

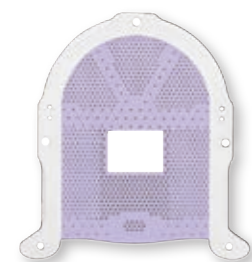

S-Shaped Head Mask
Reinforced, Openface

M24RB-VSA641W	2.4mm
M30RB-VSA671W	3.0mm


Meicen S-type Masks fit all standard S-type Baseplates from Meicen and other manufacturers.

S-type Violet Thermoplastic Masks Openface




S-Shaped Head-Neck-Shoulder Mask
Reinforced, Openface

M24RB-VSB642W	2.4mm
M30RB-VSB670W	3.0mm

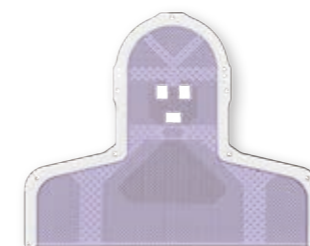


S-Shaped Head-Neck-Shoulder Mask
Reinforced, Open eyes & mouth

M24RB-VSB643W	2.4mm
M30RB-VSB672W	3.0mm

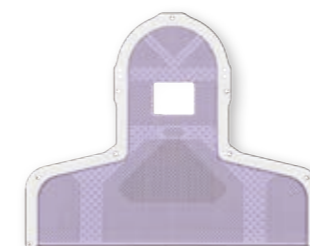

S-Shaped Head-Neck-Shoulder Mask
Reinforced, Openface

M24RB-VSB641W	2.4mm
M30RB-VSB671W	3.0mm

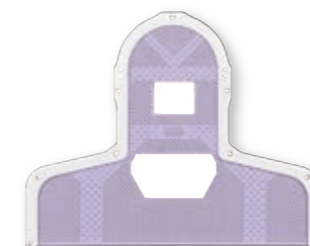

S-Shaped Head-Neck-Shoulder Mask
Reinforced, Open face & neck

M24RB-VSB644W	2.4mm
M30RB-VSB673W	3.0mm

2.0mm & 3.2mm thickness are available also, please specify thickness when requesting a quote.

L-type Violet Thermoplastic Masks



Triangular Head Mask

Standard, with nose hole

M2422-VLA191W	22% perforation, 2.4mm
M3038-VLA480W	38% perforation, 3.0mm
M2438-VLA320W	38% perforation, 2.4mm, for pediatric



Triangular Head Mask

Standard, Extended 3cm, with nose hole

M2422-VLA192W	22% perforation, 2.4mm
M3038-VLA481W	38% perforation, 3.0mm
M2415-VLA151W	15% perforation, 2.4mm



Triangular Head Mask

Reinforced, with nose hole

M24RA-VLA580W	2.4mm
M30RA-VLA610W	3.0mm



Neck-Shoulder Mask

Standard, with nose hole

M2422-VLE190W	22% perforation, 2.4mm
M3038-VLE480W	38% perforation, 3.0mm
M2422-VLE191W	22% perforation, 2.4mm, for pediatric




Neck-Thorax Mask

Standard

M2422-VLF192W	22% perforation, 2.4mm
M3038-VLF480W	38% perforation, 3.0mm


L-type Violet Thermoplastic Masks

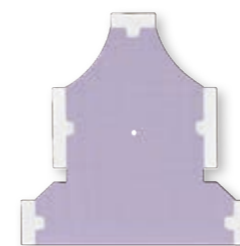

Head-Neck-Shoulder Mask

Standard, with nose hole

M2422-VLB191W	22% perforation, 2.4mm
M3038-VLB481W	38% perforation, 3.0mm
M2422-VLB192W	22% perforation, 2.4mm, for pediatric

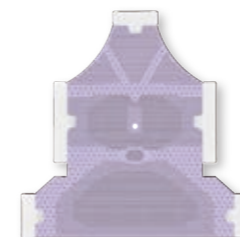

Head-Neck-Shoulder Mask

Reinforced, with nose hole


M24RA-VLB580W	2.4mm
M30RA-VLB610W	3.0mm



Chest-Pelvic Mask

Standard


M3038-VLK483W	38% perforation, 3.0mm, 228 clamp
M3038-VLK482W	38% perforation, 3.0mm, 243 clamp
M3038-VLK481W	38% perforation, 3.0mm, 247 clamp

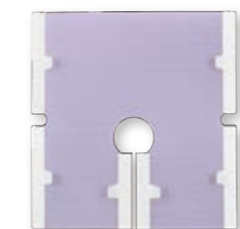

Pelvic Mask

Standard

M3038-VLL483W	38% perforation, 3.0mm, 228 clamp
M3038-VLL482W	38% perforation, 3.0mm, 243 clamp
M3038-VLL481W	38% perforation, 3.0mm, 247 clamp

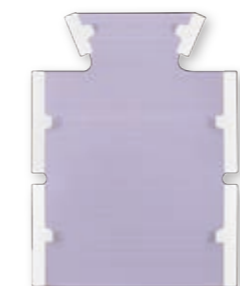

Thorax-Abdomen Mask

Standard

M2422-VLG190W	22% perforation, 2.4mm
M3038-VLG480W	38% perforation, 3.0mm

Meicen L-type Masks fit all standard O-type Baseplates from Meicen and other manufacturers.

2.0mm & 3.2mm thickness are available also, please specify thickness when requesting a quote.

L-type Violet Thermoplastic Masks Openface

P-type Violet Thermoplastic Masks



Triangular Head Mask

Reinforced, Openface


M24RB-VLA641W	2.4mm
M30RB-VLA671W	3.0mm



Triangular Head Mask

Reinforced, Open eyes & mouth

M24RB-VLA643W	2.4mm
M30RB-VLA673W	3.0mm

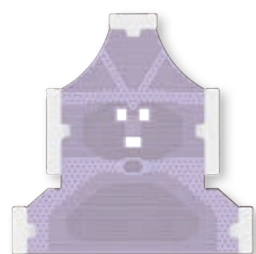

Head-Neck-Shoulder Mask

Reinforced, Open eyes & mouth


M24RB-VLB642W	2.4mm
M30RB-VLB672W	3.0mm



Head-Neck-Shoulder Mask

Reinforced, Openface

M24RB-VLB640W	2.4mm
M30RB-VLB670W	3.0mm

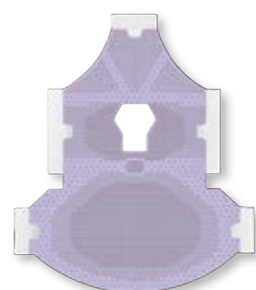

Head-Neck-Shoulder Mask

Reinforced, Extended, Openface

M24RB-VLC641W	2.4mm
M30RB-VLC671W	3.0mm



Triangular Head Mask

Standard, with nose hole

M2422-VPA190W	22% perforation, 2.4mm
M3038-VPA480W	38% perforation, 3.0mm
M2415-VPA150W	15% perforation, 2.4mm



Triangular Head Mask

Reinforced, with nose hole


M24RA-VPA580W	2.4mm
M30RA-VPA610W	3.0mm



Head-Neck-Shoulder Mask

Standard, with nose hole

M2422-VPB190W	22% perforation, 2.4mm
M3038-VPB480W	38% perforation, 3.0mm
M2415-VPB150W	15% perforation, 2.4mm

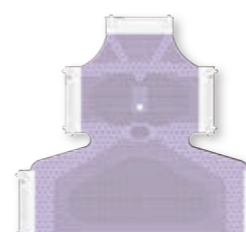

Head-Neck-Shoulder Mask

Reinforced, with nose hole

M24RA-VPB580W	2.4mm
M30RA-VPB610W	3.0mm



Chest-Pelvic Mask

Standard

M2422-VPK190W	22% perforation, 2.4mm
M3038-VPK481W	38% perforation, 3.0mm

Meicen L-type Masks fit all standard O-type Baseplates from Meicen and other manufacturers.

P-type Violet Thermoplastic Masks Openface



Triangular Head Mask

Standard, Openface

M24RB-VPA640W	2.4mm
M30RB-VPA670W	3.0mm



Triangular Head Mask

Reinforced, Openface

M24RB-VPA641W	2.4mm
M30RB-VPA671W	3.0mm



Triangular Head Mask

Reinforced, Open eyes & mouth

M24RB-VPA642W	2.4mm
M30RB-VPA672W	3.0mm



Head-Neck-Shoulder Mask

Reinforced, Openface

M24RB-VPB640W	2.4mm
M30RB-VPB670W	3.0mm



Head-Neck-Shoulder Mask

Reinforced, Open eyes & mouth

M24RB-VPB641W	2.4mm
M30RB-VPB671W	3.0mm

Pin-lock type Violet Thermoplastic Masks

Compatible with Meicen C & Y-Series Baseplates



Neck-Thorax Mask

Reinforced

M24RA-VSF580W	2.4mm
M30RA-VSF610W	3.0mm



Trunk Mask (Chest-Pelvic)

Standard

M2422-VST190W	22% perforation, 2.4mm
M3038-VST480W	38% perforation, 3.0mm



Trunk Mask (Pelvic)

Standard, with Groin Lock

M2422-VST191W	22% perforation, 2.4mm
M3038-VST481W	38% perforation, 3.0mm

2.0mm & 3.2mm thickness are available also, please specify thickness when requesting a quote.

H-type Violet Thermoplastic Masks



Trunk Mask (Chest-Pelvic)

Standard

M2422-VHT190W	22% perforation, 2.4mm
M3038-VHT480W	38% perforation, 3.0mm



Trunk Mask (Pelvic)

Standard, with Groin Lock

M2422-VHT191W	22% perforation, 2.4mm
M3038-VHT481W	38% perforation, 3.0mm



Trunk Mask (Chest-Pelvic)

Reinforced

M24RA-VHT580W	2.4mm
M30RA-VHT610W	3.0mm



Trunk Mask (Chest-Pelvic)


Solid

M2400-VHT120W	2.4mm
M3000-VHT420W	3.0mm


2.0mm & 3.2mm thickness are available also, please specify thickness when requesting a quote.

X-Knife Violet Thermoplastic Masks

Meicen X-Knife Violet Thermoplastic Masks are composed of a new, premium thermoplastic material designed specifically to offer the ultimate in patient comfort with exceptional rigidity and extremely low shrinkage. Additionally, the X-Knife Violet Thermoplastic Masks have a non-stick coating applied to the patient-facing mask surfaces, helping ensure that the masks do not stick to the patient's hair, beard, etc.

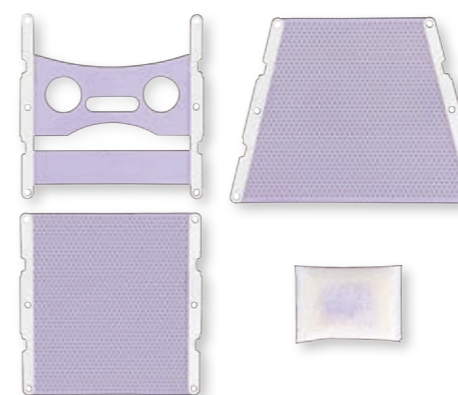


X-Knife Mask
M3038-VXA481W

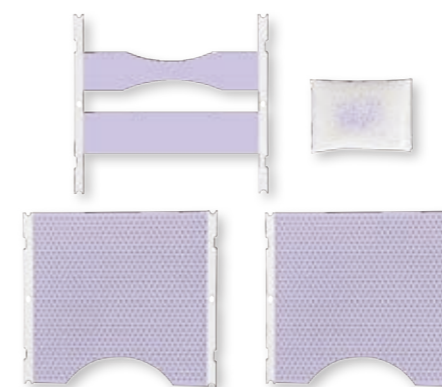

X-Knife Mask
M3038-VXA490W

- Easy to manipulate and mould
- Adhesive Layers
- Low shrinkage
- Cost-effective

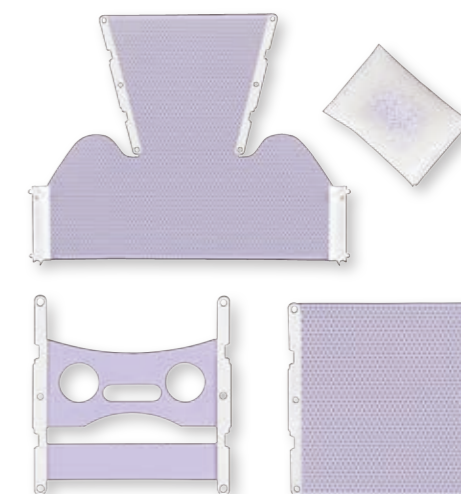

X-Knife Head-Shoulder Mask
M3038-VXB480W

E-type Violet Thermoplastic Masks



Head Mask

Standard

M2422-VEA190W	22% perforation, 2.4mm
M3038-VEA480W	38% perforation, 3.0mm
M2415-VEA150W	15% perforation, 2.4mm



Head Mask

Reinforced

M24RA-VEA580W	2.4mm
M30RA-VEA610W	3.0mm



Head Mask

Standard, Extended 5cm


M2422-VEA192W	22% perforation, 2.4mm
M3038-VEA482W	38% perforation, 3.0mm



Head Mask

Standard, Extended 21cm

M2422-VEA193W	22% perforation, 2.4mm
M3038-VEA483W	38% perforation, 3.0mm

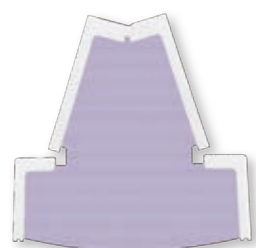


Head-Neck-Shoulder Mask

Standard

M2422-VEB190W	22% perforation, 2.4mm
M3038-VEB480W	38% perforation, 3.0mm
M2415-VEB150W	15% perforation, 2.4mm

E-type Violet Thermoplastic Masks

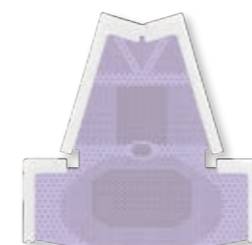

Head-Neck-Shoulder Mask

Reinforced

M24RA-VEB580W	2.4mm
M30RA-VEB610W	3.0mm



Chest-Pelvic Mask

Standard

M3038-VEK480W	38% perforation, 3.0mm
---------------	------------------------



Triangular Head Mask

Reinforced

M24RA-VEA600W	2.4mm
M30RA-VEA620W	3.0mm



Triangular Head Mask

Reinforced

M24RA-VEA602W	2.4mm
M30RA-VEA621W	3.0mm

2.0mm & 3.2mm thickness are available also, please specify thickness when requesting a quote.

V-type Violet Thermoplastic Masks



V-Shaped Head Mask

Standard

M2422-VVA190W	22% perforation, 2.4mm
M3038-VVA480W	38% perforation, 3.0mm



V-Shaped Head Mask

Reinforced

M24RA-VVA580W	2.4mm
M30RA-VVA610W	3.0mm



V-Shaped Head Mask

Reinforced, Open eyes & mouth

M24RB-VVA640W	2.4mm
M30RB-VVA670W	3.0mm



V-Shaped Head Mask

Reinforced, Openface

M24RB-VVA641W	2.4mm
M30RB-VVA671W	3.0mm




V-Shaped Head Mask

Standard, Extended

Extended 5cm	Extended 9cm	
M2422-VVA191W	M2422-VVA192W	22% perforation, 2.4mm
M3038-VVA481W	M3038-VVA482W	38% perforation, 3.0mm


V-type Violet Thermoplastic Masks

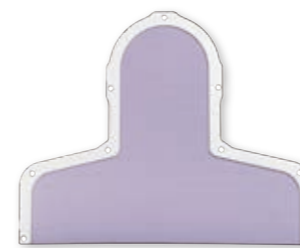

V-Shaped Head-Neck-Shoulder Mask

Standard

M2422-VVB190W	22% perforation, 2.4mm
M3038-VVB480W	38% perforation, 3.0mm

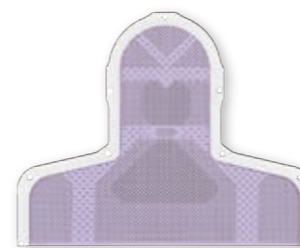

V-Shaped Head-Neck-Shoulder Mask

Reinforced


M24RA-VVB580W	2.4mm
M30RA-VVB610W	3.0mm



V-Shaped Head-Neck-Shoulder Mask

Reinforced, Openface

M24RB-VVB640W	2.4mm
M30RB-VVB670W	3.0mm

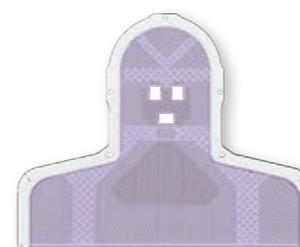


V-Shaped Head-Neck-Shoulder Mask

Reinforced, Open eyes & mouth

M24RB-VVB641W	2.4mm
M30RB-VVB671W	3.0mm

Meicen Violet Thermoplastic Sheets

Meicen Violet Thermoplastic Sheets offer flexible solution for radiation therapy immobilization, available in a wide variety of thickness, perforation and dimension.

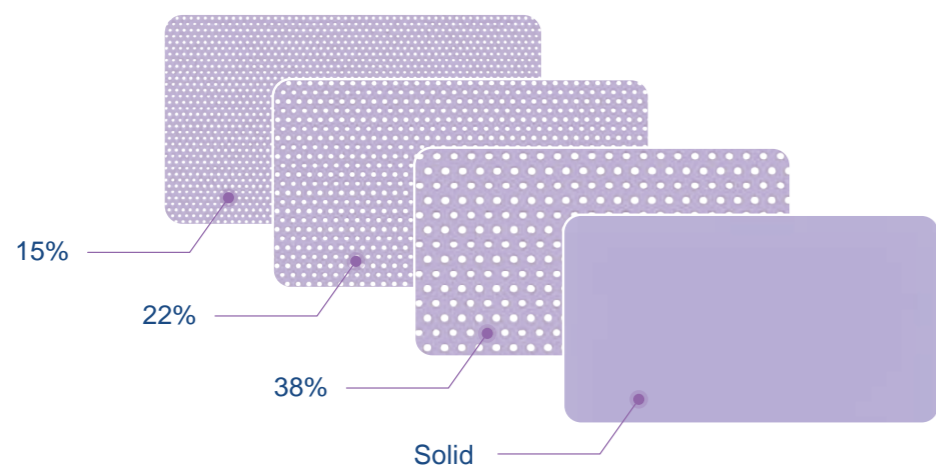

REF No.	Dimension	Thickness	Perforation
M2415-VRT150	900x600x2.4mm	2.4mm	15%
M2422-VRT190	900x600x2.4mm	2.4mm	22%
M2438-VRT320	900x600x2.4mm	2.4mm	38%
M3038-VRT480	900x580x3.0mm	3.0mm	38%
M2422-VRT191	600x450x2.4mm	2.4mm	22%
M3038-VRT481	600x450x3.0mm	3.0mm	38%
M2400-VRT121	600x450x2.4mm	2.4mm	Solid
M3000-VRT420	600x450x3.0mm	3.0mm	Solid

2.0mm & 3.2mm thickness are available also, please specify thickness when requesting a quote.

Meicen Diagnostic Couchtop Overlays

- The Couchtop Overlays provide a clean, flat surface for accurate and easily reproducible patient positioning.
- The internal density distribution of the product is uniform, and there is no unacceptable artifact in the scanning area.
- Can be quickly installed and disassembled, simple and convenient.
- It is stable and reliable after installation and meets the leveling requirements.
- The Couchtop Overlays are designed to work with Indexing bars as well as other positioning systems.

* The Couchtop Overlays can be customized according to the user's CT or MRI equipment.

MC-CTF019 CT Overlay



Also Available As MR Compatible

MC-MRG002

MR

